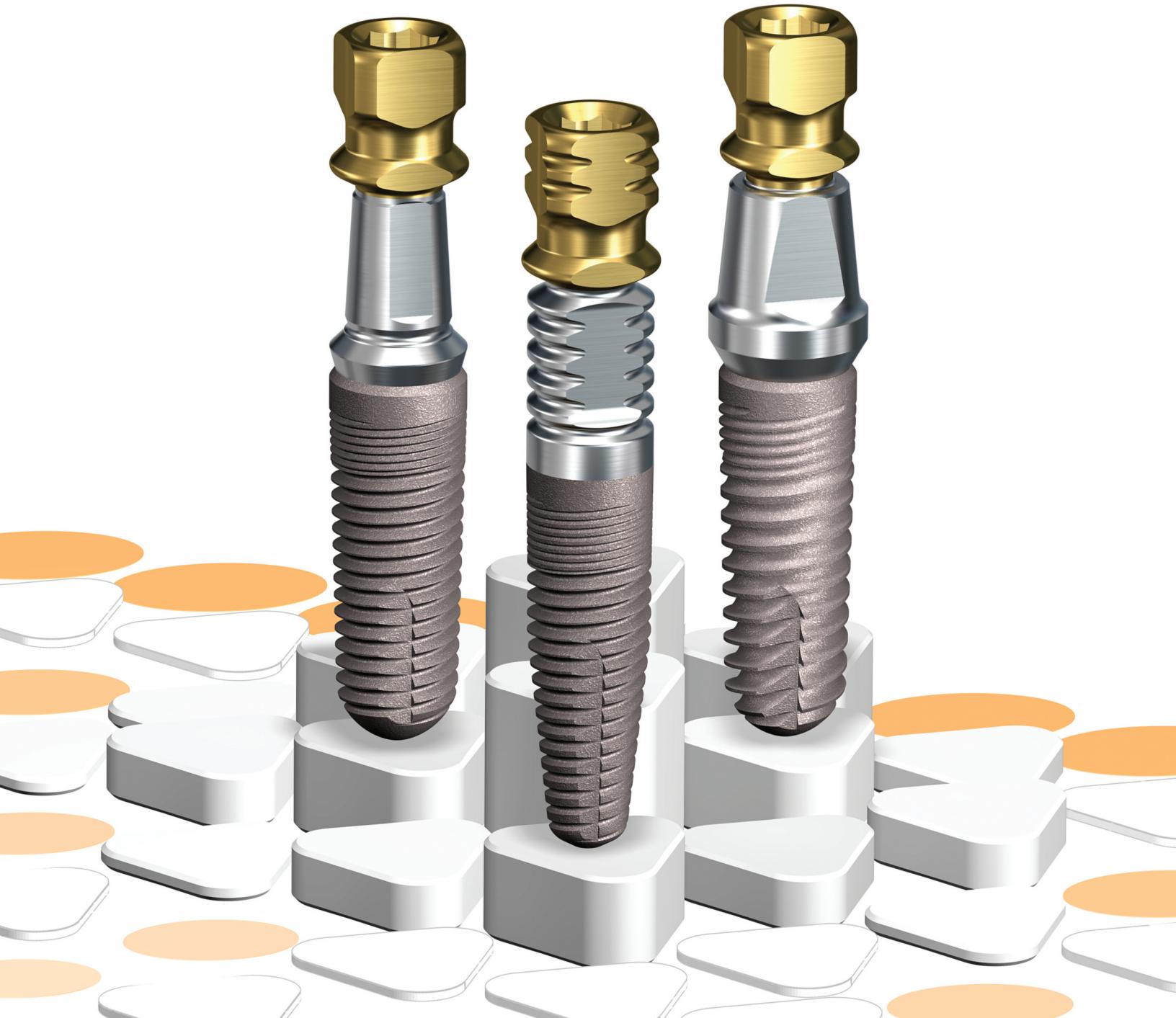





RePlant® System Catalog





Mini-threads with Double-lead Threads for Faster Insertion

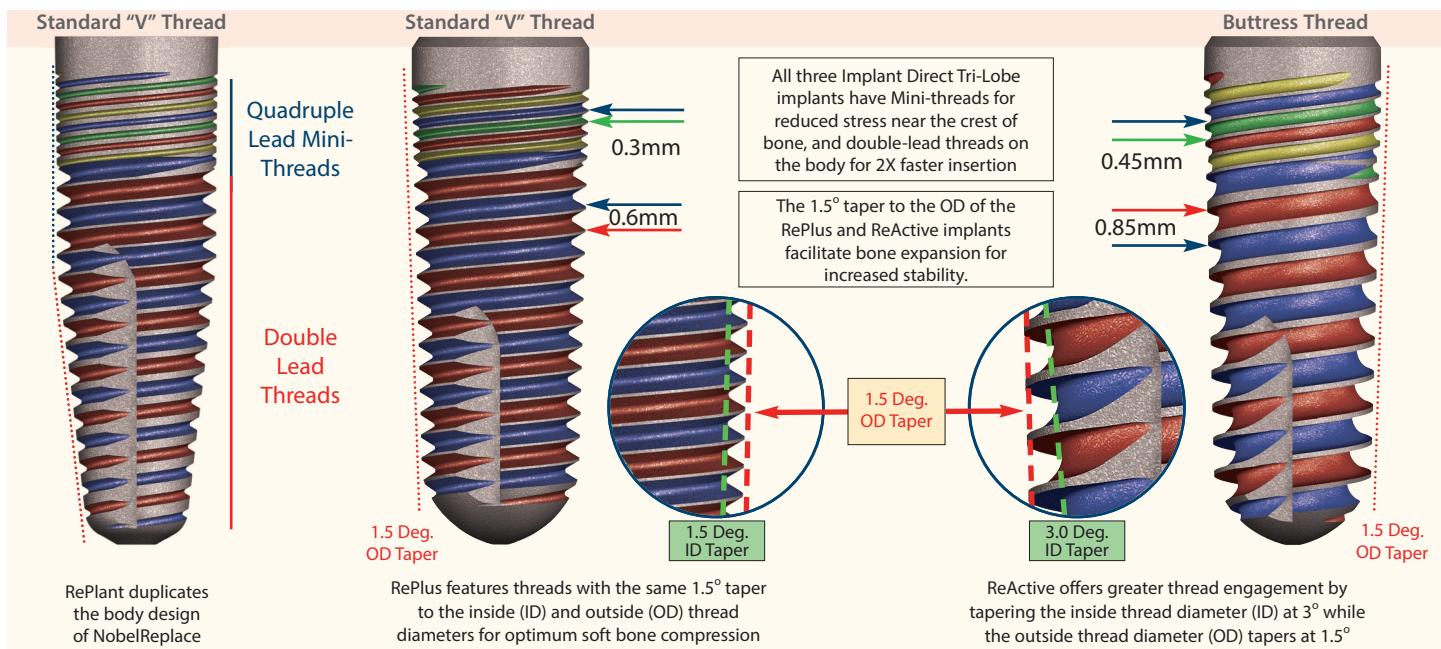The bodies of the RePlus and ReActive implants are tapered evenly from the top to facilitate insertion into an undersized socket in soft bone providing compression that increases initial stability (G. Niznick: *Clinical Oral Health 2000 "Achieving Success in Soft Bone"*). The RePlant implant is designed with a straight/tapered profile matching that of the NobelReplace Implant (Nobel Biocare). Mini-threads start 1.0mm below the crest of the bone and extend approximately 2.0mm before transitioning into double-lead threads to the apex of the implant. The mini-threads, half the depth and pitch of the double-lead body threads, reduce stress in the critical crestal bone region, while double-lead threads reduce the number of turns required to seat the implant. Vertical cutting grooves extend halfway up the implant from the apex for self-tapping insertion. All three implants are available with SBM textured surface over the entire implant and HA is an option for RePlus and ReActive.




SBM™ and HA Coated Surface Options


SBM surface is created by using a **Soluble Blast Media** of Hydroxylapatite (HA) crystals to create a medium-rough texture. We use the same vendor for its SBM and HA coated surface treatments as used by Paragon Implant Company from 1992-2000 when these studies were conducted.

SEM's of Our Surfaces

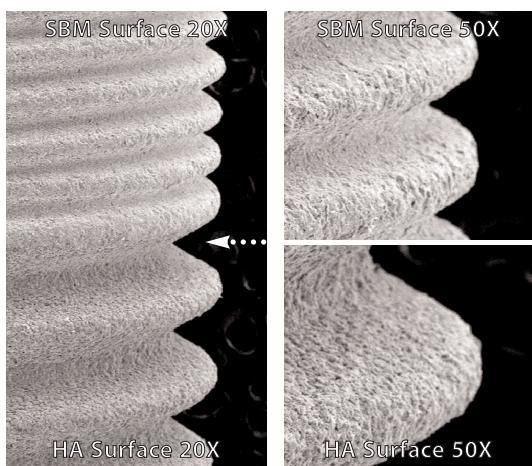
Mini-threads textured by blasting with SBM™: Blasted with Soluble HA Particles
HA Coating: High Crystalline 40u layer
RePlus & ReActive Implants have an HA option.



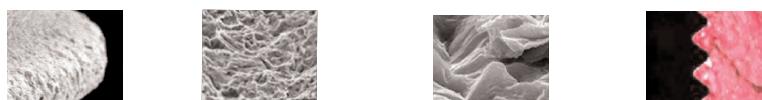Animal Studies Bone Response to SBM Surface



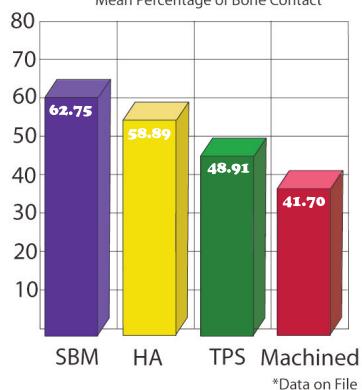Brazilian Study* Shows SBM Provides Superior Bone Contact
Mean Percentage of Bone Contact



LSU Study* Shows SBM Provides Superior Bone Contact
Mean Percentage of Bone Contact

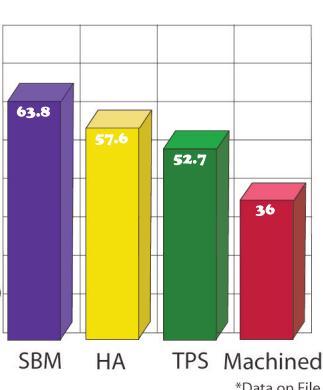


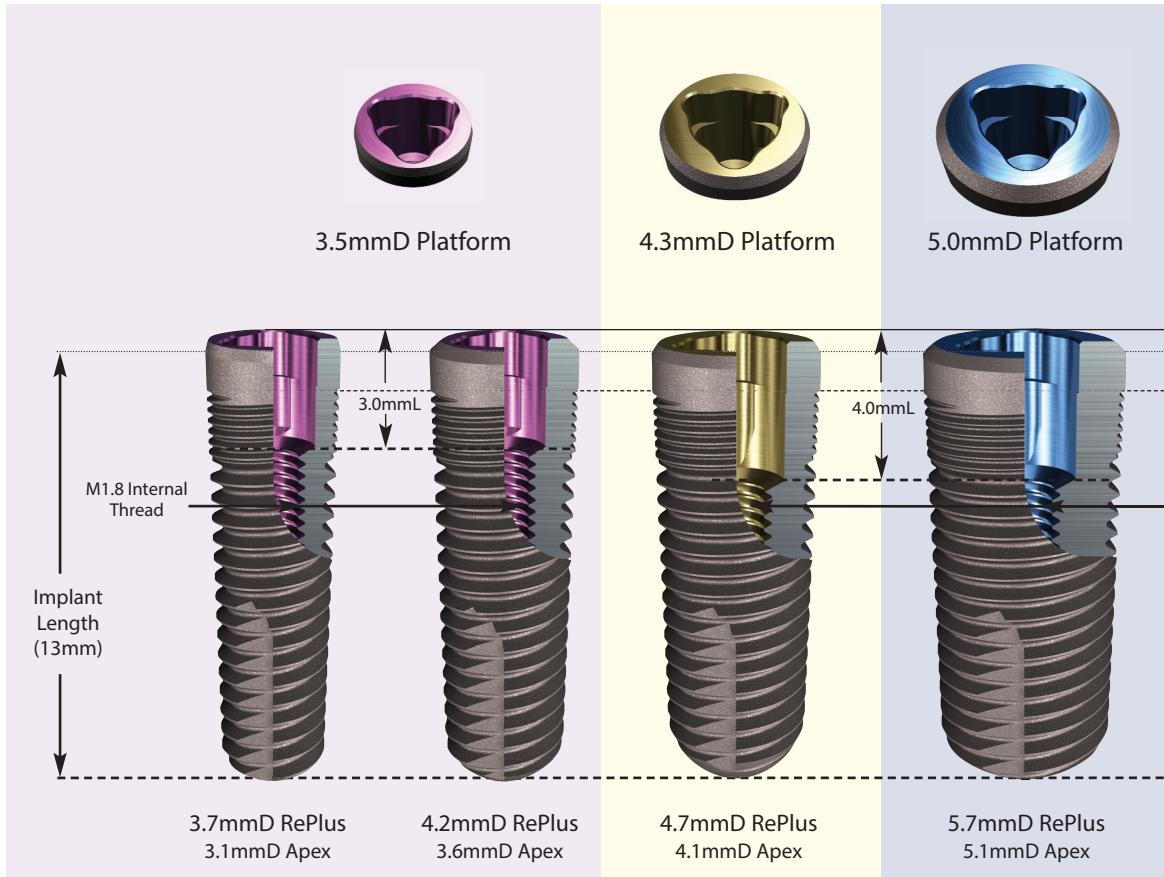


Table of Contents


	Page #	
	Implant Design & Surface	1
		2-3
	RePlant, RePlus & ReActive Implants	4-9
		10
	Healing Collars, Transfers, Analogs & Alignment Tools	11
		12
	Straight Abutment & Ancillary Components	13
		14
	Abutments for Partially Edentulous Restorations	15-17
		18-19
	Procedures for Titanium FMT / Temporary Abutment	20
		21-24
	Surgical & Prosthetic Instruments	25
		26-29
	RePlant Surgical Sequence & Drills	30-33
Implant Dimensions & Product Label Description		
		
	Selection Chart for Implant & Straight Abutment Components	
		
	Procedures for Contoured Abutments & Lab Alignment Tools	
	Screw-Receiving, Ball, GPS Abutments & Components	
	RePlus & ReActive Surgical Sequence & Drills	




RePlus and ReActive Implants




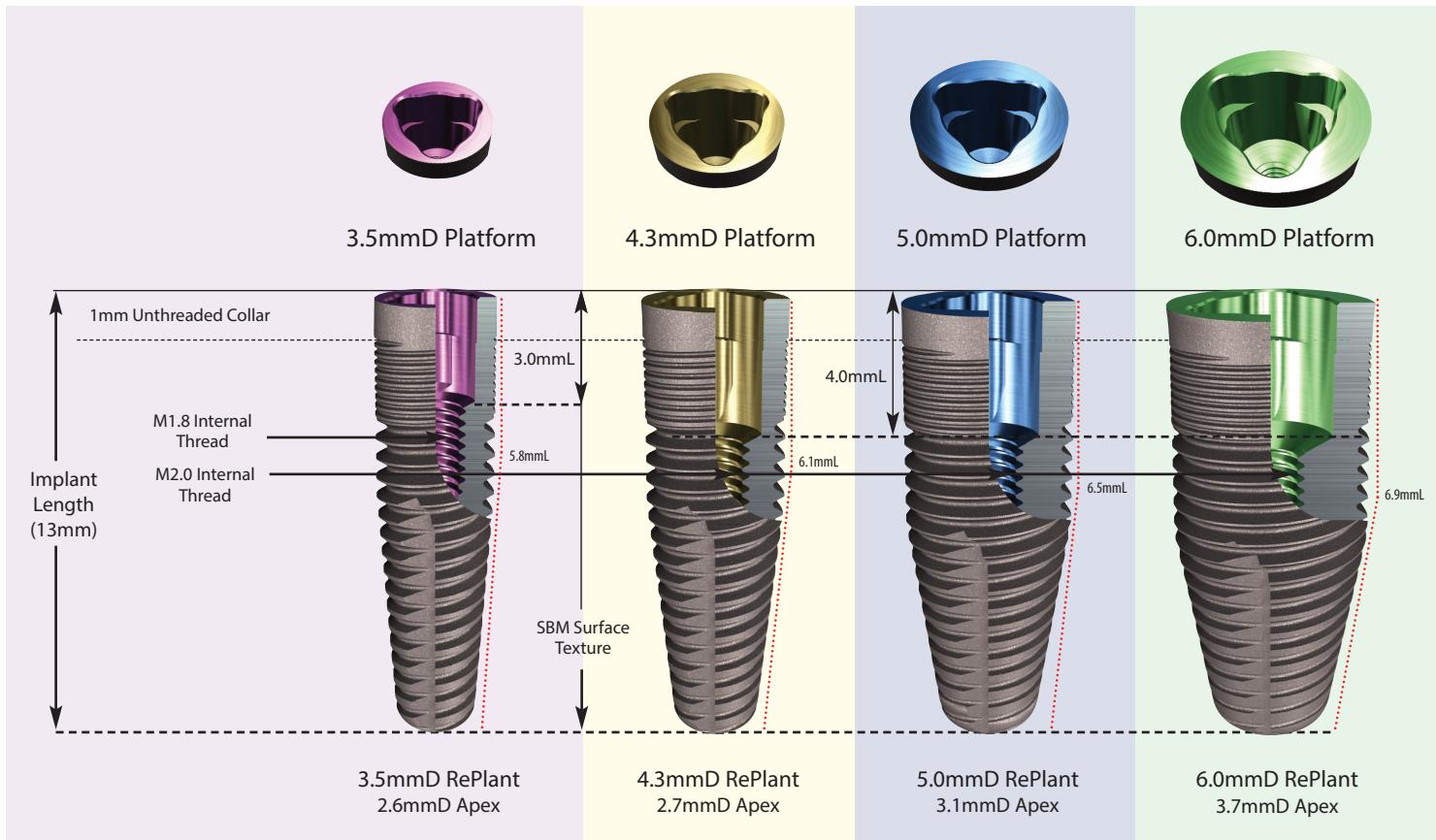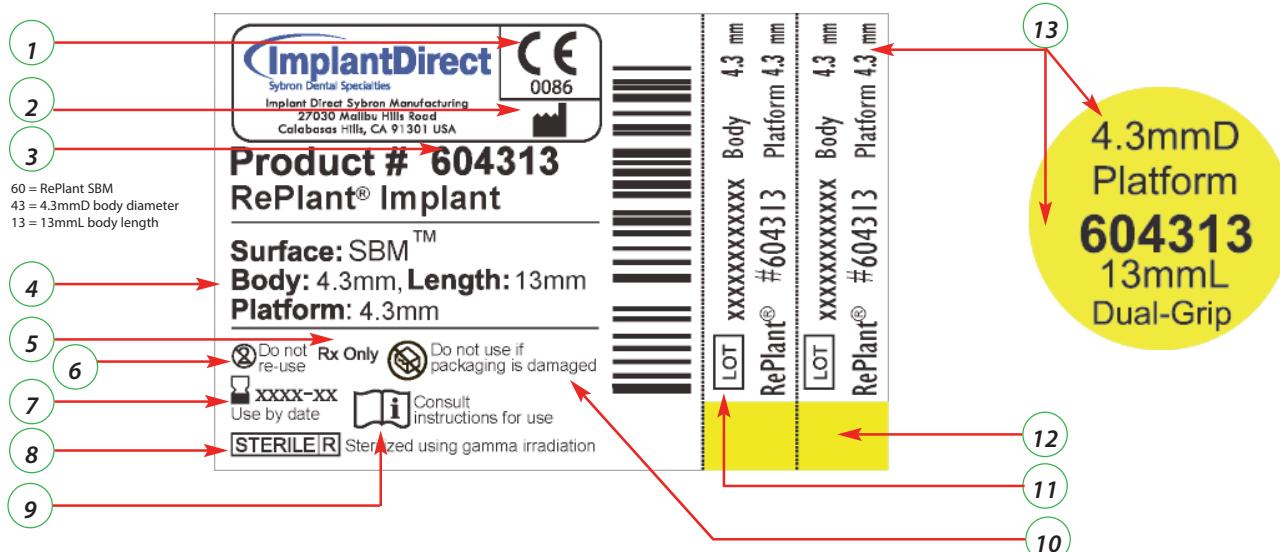RePlant Implants



RePlant *Implant* Vial and Cap Label

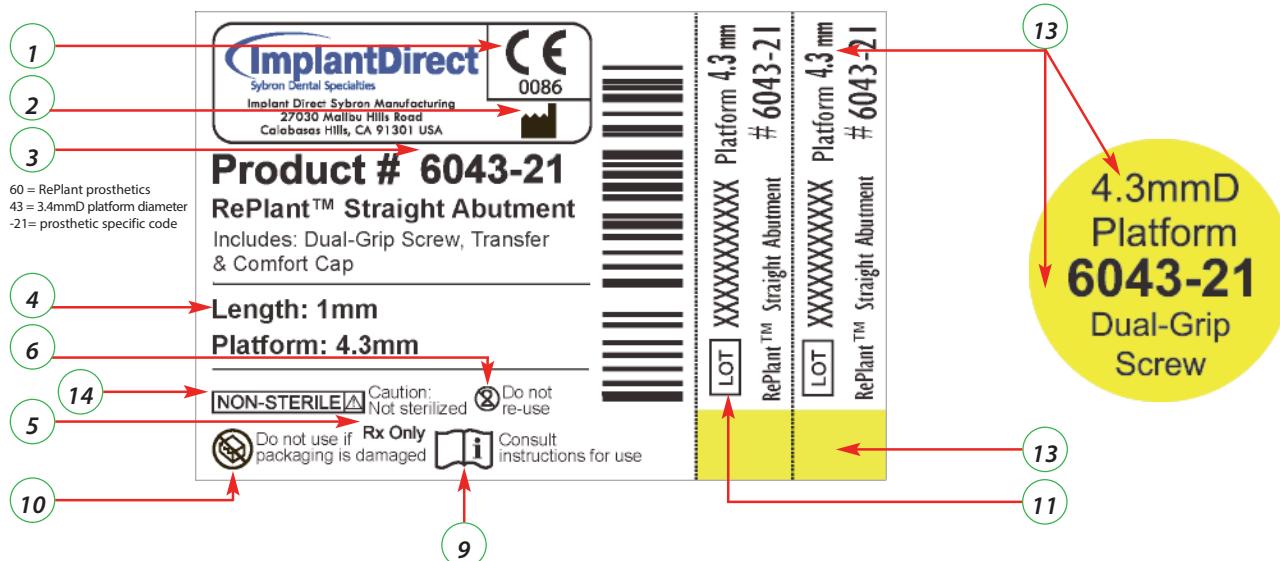RePlant Prosthetic Vial and Cap Label




- | | |
|-------------------------|--|
| CE | 1) Our products carry the CE mark and fulfill the requirements of the Medical Device Directive 93/42/EEC |
| [] | 2) Manufacturer Symbol |
| | 3) Product code or reference number |
| | 4) Product description |
| Rx Only | 5) Caution: US law restricts the sale of these devices by or on the order of a dentist or physician |
| (X) | 6) Do not re-use symbol |
| XXXX-XX | 7) Use by date (XXXX-XX) |
| STERILE[R] | 8) Sterilized using gamma irradiation symbol |
| [i] | 9) Refer to Instructions For Use (IFU) |
| [] | 10) Tamper proof symbol |
| | 11) Lot number (XXXXXXXXXX) |
| | 12) Implant body diameter color reference on surgical tray |
| | 13) Implant prosthetic platform diameter (written and color-coded) |
| [NON-STERILE[R]] | 14) Non-sterile symbol |



RePlant



RePlant Two-Piece Implant System Featuring:

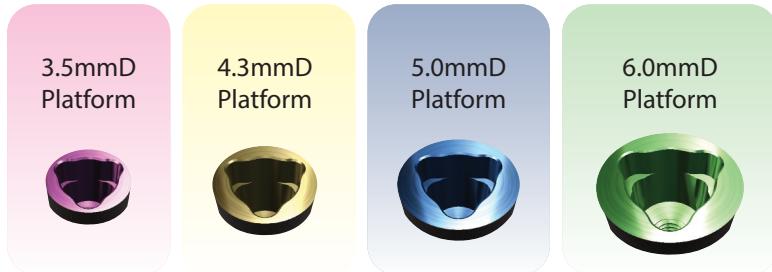- Straight/Tapered body design and dimensions match NobelReplace® Tapered Groovy implant allowing insertion with Nobel Biocare's drills
- Tri-Lobe color-coded platform is compatible with NobelReplace® abutments
- SBM (Soluble Blast Media) textured surface to top of implant
- Color-coded fixture-mount usable as a transfer and shortened for use as a temporary abutment
- Use of titanium alloy for an increase in strength compared to Nobel Biocare's Tri-Lobe implants manufactured from CP titanium



Carrier/Transfer

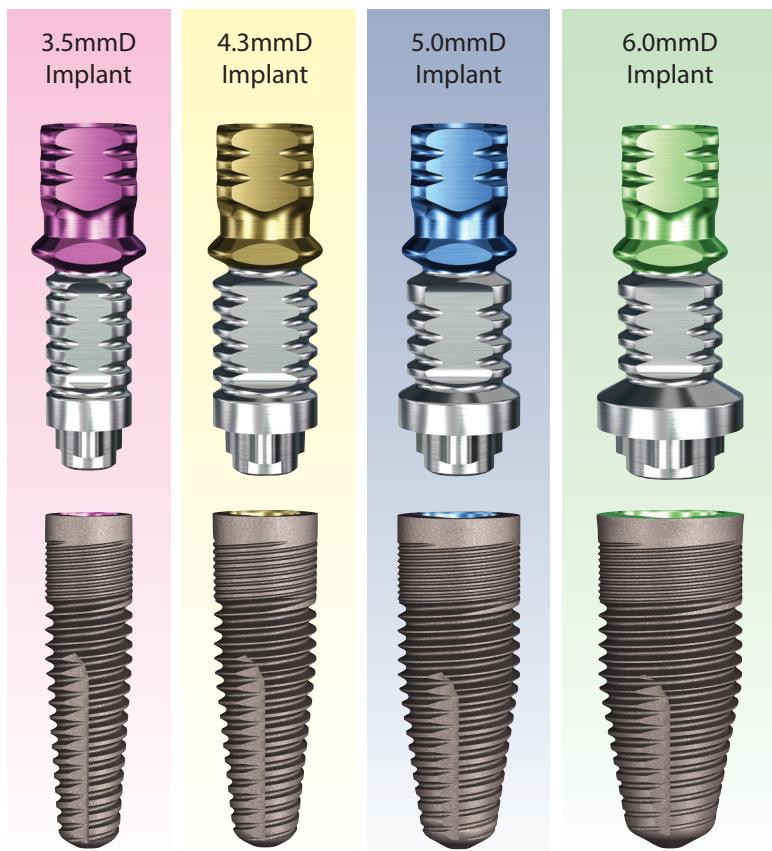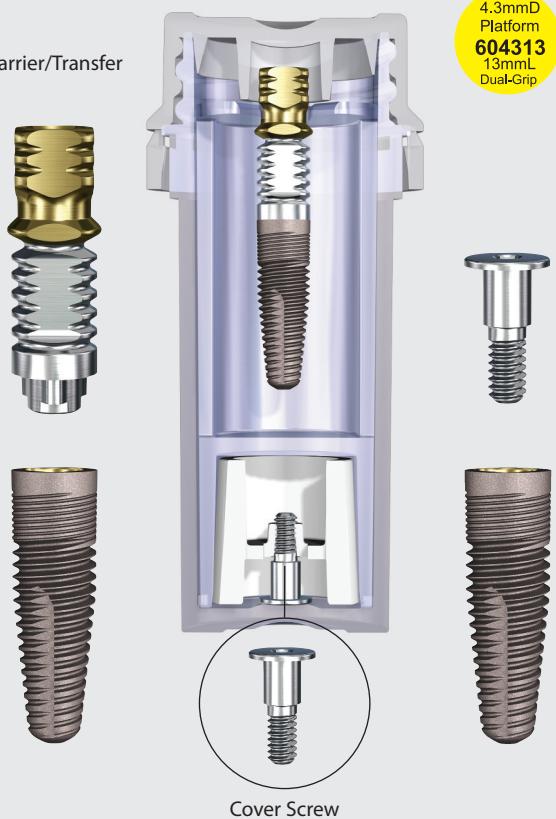4.3mmD
Platform
604313
13mmL
Dual-Grip



Placing Carrier/Transfers & Making Impressions



Color-coded Titanium
fixture-mount carrier



Align flat to labial for
indexing of tri-lobe



Fixture-mount is
screw-retained



Use fixture-mount as
closed tray transfer



RePlant


RePlant: Two-Piece Implant

**USE REPLANT SURGICAL TRAY OR NOBEL BIOCARE'S TAPERED
REPLACE IMPLANT DRILLS FOR INSERTION OF REPLANT IMPLANT**


SBM Surface Implants

3.5mmD*Includes Fixture Mount, Dual-Grip™ Fixation and Cover Screw*

Code	Diameter	Length	Platform	Price
603508	3.5mmD	8mmL	3.5mmD	\$175.00
603510	3.5mmD	10mmL	3.5mmD	\$175.00
603511	3.5mmD	11.5mmL	3.5mmD	\$175.00
603513	3.5mmD	13mmL	3.5mmD	\$175.00
603516	3.5mmD	16mmL	3.5mmD	\$175.00

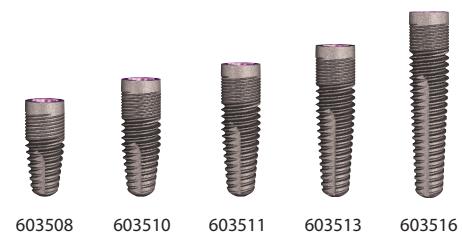**4.3mmD**

Code	Diameter	Length	Platform	Price
604308	4.3mmD	8mmL	4.3mmD	\$175.00
604310	4.3mmD	10mmL	4.3mmD	\$175.00
604311	4.3mmD	11.5mmL	4.3mmD	\$175.00
604313	4.3mmD	13mmL	4.3mmD	\$175.00
604316	4.3mmD	16mmL	4.3mmD	\$175.00

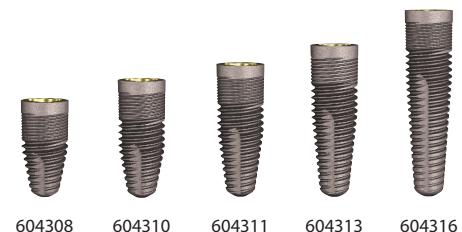**5.0mmD**

Code	Diameter	Length	Platform	Price
605008	5.0mmD	8mmL	5.0mmD	\$175.00
605010	5.0mmD	10mmL	5.0mmD	\$175.00
605011	5.0mmD	11.5mmL	5.0mmD	\$175.00
605013	5.0mmD	13mmL	5.0mmD	\$175.00
605016	5.0mmD	16mmL	5.0mmD	\$175.00

**6.0mmD**

Code	Diameter	Length	Platform	Price
606008	6.0mmD	8mmL	6.0mmD	\$175.00
606010	6.0mmD	10mmL	6.0mmD	\$175.00
606011	6.0mmD	11.5mmL	6.0mmD	\$175.00
606013	6.0mmD	13mmL	6.0mmD	\$175.00
606016	6.0mmD	16mmL	6.0mmD	\$175.00



Replacement Surgical Cover Screws

*Components include Dual-Grip™ feature for use with
1.25mmD Hex Tool or Nobel Unigrip ScrewDriver*

Code	Platform	Price
6035-01	3.5mmD	\$20.00
6043-01	4.3mmD	\$20.00
6050-01	5.0mmD	\$20.00
6060-01	6.0mmD	\$20.00

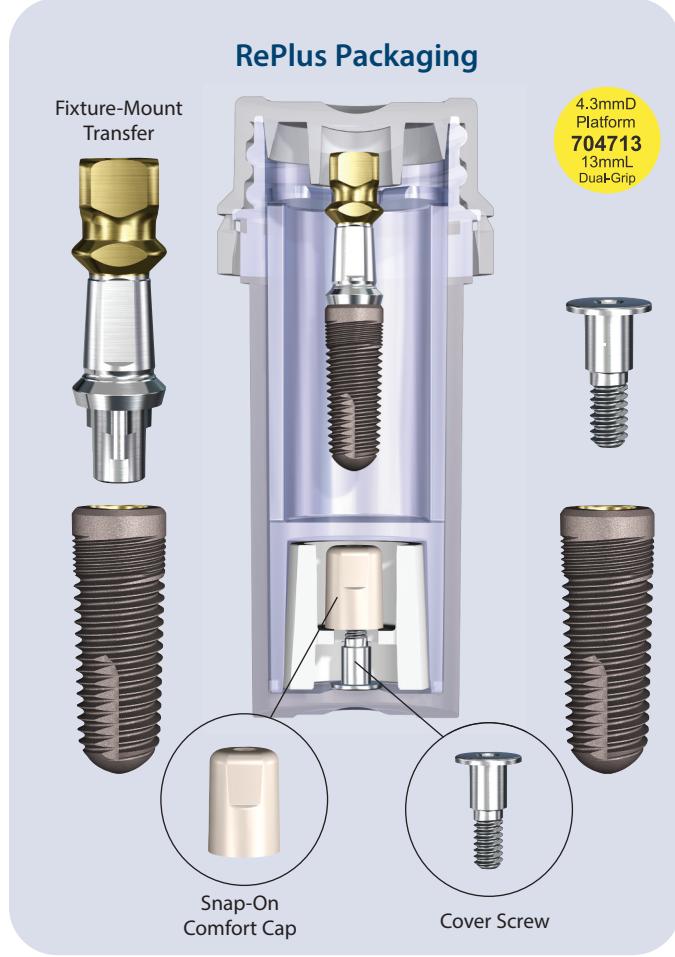




RePlus

RePlus Two-Piece Implant System Featuring:

- Tri-Lobe color-coded platform is compatible with NobelReplace® abutments
- Color-coded fixture-mount usable as a transfer and shortened for use as a straight final abutment with snap-on comfort cap
- Surface options: SBM (Soluble Blast Media) or HA-coated
- Platform-switching with RePlant and NobelReplace abutments



Fixture-Mount Transfer: Shortened for final abutment with Snap-on Comfort Cap and Transfer





RePlus: Two-Piece Implant

RePlus

SBM and HA Surface Implants Includes Fixture Mount with Dual-Grip™ Fixation Screw, Cover Screw and Comfort Cap

3.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
703708	\$200.00	713708	\$225.00	3.7mmD	8mmL	3.5mmD
703710	\$200.00	713710	\$225.00	3.7mmD	10mmL	3.5mmD
703711	\$200.00	713711	\$225.00	3.7mmD	11.5mmL	3.5mmD
703713	\$200.00	713713	\$225.00	3.7mmD	13mmL	3.5mmD
703716	\$200.00	713716	\$225.00	3.7mmD	16mmL	3.5mmD



4.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
704208	\$200.00	714208	\$225.00	4.2mmD	8mmL	3.5mmD
704210	\$200.00	714210	\$225.00	4.2mmD	10mmL	3.5mmD
704211	\$200.00	714211	\$225.00	4.2mmD	11.5mmL	3.5mmD
704213	\$200.00	714213	\$225.00	4.2mmD	13mmL	3.5mmD
704216	\$200.00	714216	\$225.00	4.2mmD	16mmL	3.5mmD



4.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
704708	\$200.00	714708	\$225.00	4.7mmD	8mmL	4.3mmD
704710	\$200.00	714710	\$225.00	4.7mmD	10mmL	4.3mmD
704711	\$200.00	714711	\$225.00	4.7mmD	11.5mmL	4.3mmD
704713	\$200.00	714713	\$225.00	4.7mmD	13mmL	4.3mmD
704716	\$200.00	714716	\$225.00	4.7mmD	16mmL	4.3mmD



5.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
705708	\$200.00	715708	\$225.00	5.7mmD	8mmL	5.0mmD
705710	\$200.00	715710	\$225.00	5.7mmD	10mmL	5.0mmD
705711	\$200.00	715711	\$225.00	5.7mmD	11.5mmL	5.0mmD
705713	\$200.00	715713	\$225.00	5.7mmD	13mmL	5.0mmD
705716	\$200.00	715716	\$225.00	5.7mmD	16mmL	5.0mmD

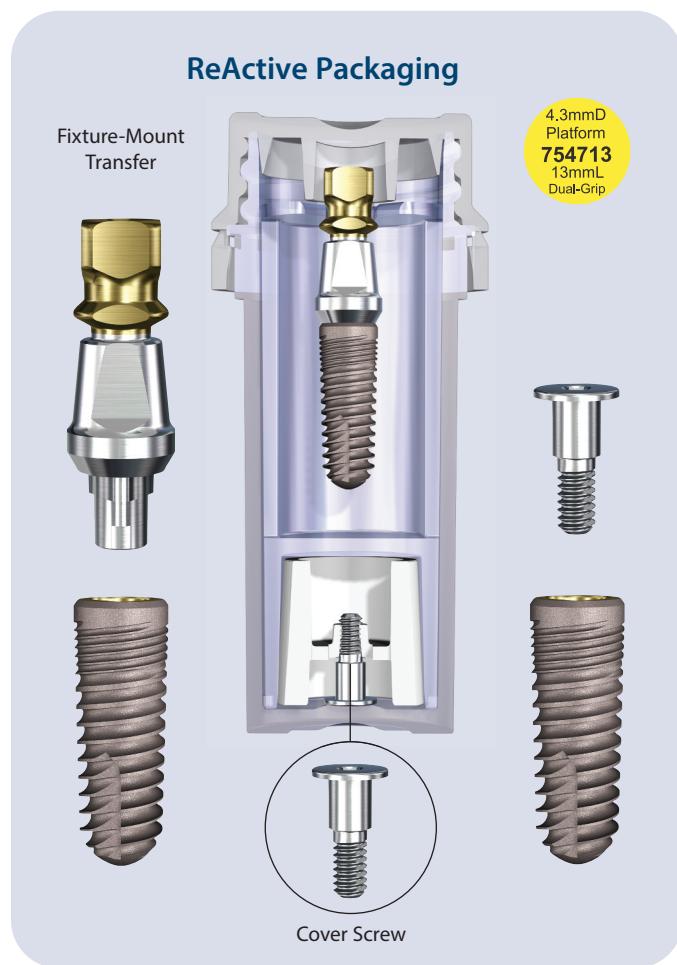
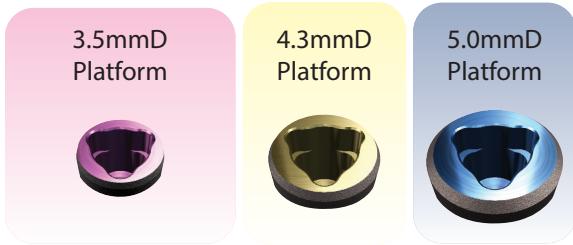





ReActive

ReActive Two-Piece Implant System Featuring:


- Tri-Lobe color-coded platform is compatible with NobelReplace® abutments
- Color-coded fixture-mount usable as a transfer and shortened for use as a platform-switching straight final abutment
- Surface options: SBM (Soluble Blast Media) or HA-coated
- Platform-switching with RePlant and NobelReplace abutments
- Buttress threads with increased thread depth providing greater surface area for bone attachment and primary implant stability



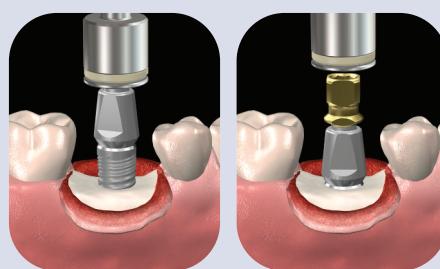Placing Carrier/Transfers & Taking Impressions



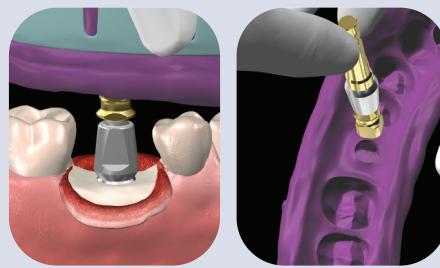
Color-coded titanium fixture-mount carrier
Align flat to labial for indexing of tri-lobe



Fixture-mount is screw retained
Use fixture-mount as closed tray transfer

Fixture-Mount Transfer:
Shortened to fabricate final abutment

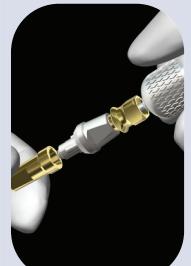Step 1



Step 2





ReActive: Two-Piece Implant

ReActive

SBM and HA Surface Implants Includes Fixture Mount with Dual-Grip™ Fixation Screw and Cover Screw

3.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
753708	\$200.00	763708	\$225.00	3.7mmD	8mmL	3.5mmD
753710	\$200.00	763710	\$225.00	3.7mmD	10mmL	3.5mmD
753711	\$200.00	763711	\$225.00	3.7mmD	11.5mmL	3.5mmD
753713	\$200.00	763713	\$225.00	3.7mmD	13mmL	3.5mmD
753716	\$200.00	763716	\$225.00	3.7mmD	16mmL	3.5mmD



4.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
754208	\$200.00	764208	\$225.00	4.2mmD	8mmL	3.5mmD
754210	\$200.00	764210	\$225.00	4.2mmD	10mmL	3.5mmD
754211	\$200.00	764211	\$225.00	4.2mmD	11.5mmL	3.5mmD
754213	\$200.00	764213	\$225.00	4.2mmD	13mmL	3.5mmD
754216	\$200.00	764216	\$225.00	4.2mmD	16mmL	3.5mmD



4.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
754708	\$200.00	764708	\$225.00	4.7mmD	8mmL	4.3mmD
754710	\$200.00	764710	\$225.00	4.7mmD	10mmL	4.3mmD
754711	\$200.00	764711	\$225.00	4.7mmD	11.5mmL	4.3mmD
754713	\$200.00	764713	\$225.00	4.7mmD	13mmL	4.3mmD
754716	\$200.00	764716	\$225.00	4.7mmD	16mmL	4.3mmD



5.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
755708	\$200.00	765708	\$225.00	5.7mmD	8mmL	5.0mmD
755710	\$200.00	765710	\$225.00	5.7mmD	10mmL	5.0mmD
755711	\$200.00	765711	\$225.00	5.7mmD	11.5mmL	5.0mmD
755713	\$200.00	765713	\$225.00	5.7mmD	13mmL	5.0mmD
755716	\$200.00	765716	\$225.00	5.7mmD	16mmL	5.0mmD





Fixture-Mount Packaging

The RePlant, RePlus and ReActive implants are mounted on color-coded fixture-mounts, designed for use as implant level transfers or shortened to function as abutments. A cover screw is included in each vial. The RePlus fixture-mount is provided with a snap-on comfort cap.



Prepare surgical site with double cutting step drill



All the implants come in a sterile vial



Each label is color-coded with two tear-off labels



Each cap is color-coded by platform diameter of implant



Drop implant assembly out of vial onto a sterile surface



Remove implant assembly with a fixture-mount drill or ratchet



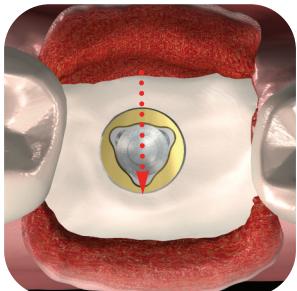
Insert implant into osteotomy



Index fixture-mount flat to labial



Remove fixture-mount



The flat of the fixture mount aligned the lobe to the labial



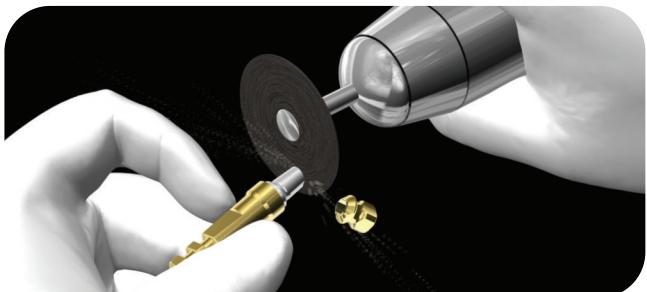
Remove the holder containing the cover screw & comfort cap



Snap-on comfort cap is located in the top portion of holder



Attach fixture-mount to corresponding color-coded implant analog



Remove colored section from fixture-mount through groove to create a final straight abutment



Attach abutment using 30Ncm in hard bone & finger pressure in soft bone



Check for occlusal clearance prior to attaching comfort cap



Test retention & rotate if necessary to increase retention



Suture soft tissue around comfort cap



The flat side of the Fixture Mount/Transfer [Fig. 1] and the angle of the 15° abutment [Fig. 2] are manufactured to align with the outer radius of one of the implant's three internal lobes.

Orienting the 15° Angled Abutment

Fig. 2 The two-piece, 15° Angled Abutment consists of an angled body and a fixation screw




To establish the optimum angle of the abutment, index the flat surface of the Fixture Mount/Transfer either toward or away from the desired direction of abutment angle at the time of implant placement.

Option 1: Attaching an Angled Abutment and chairside procedures for preparation



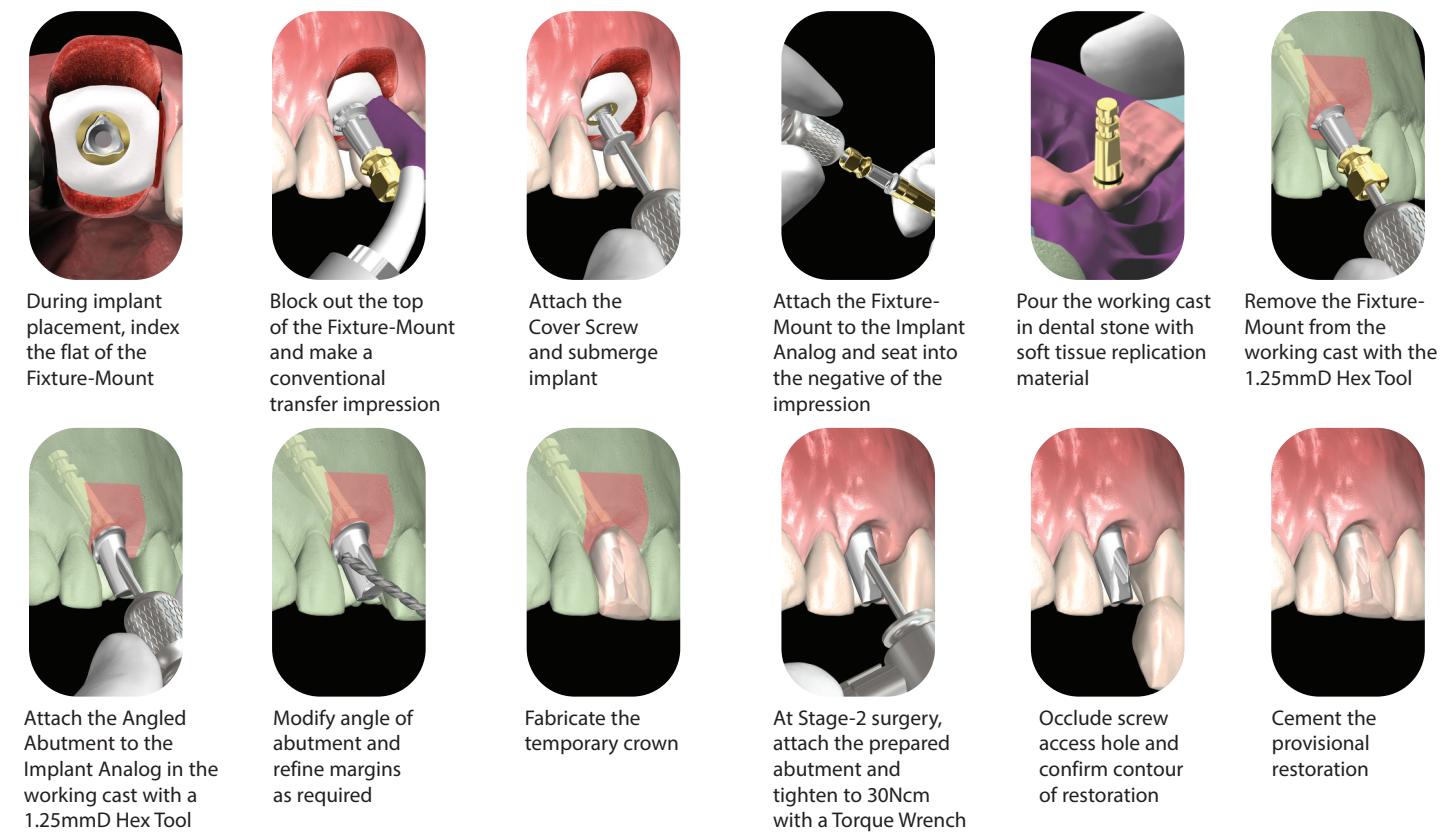Option 2: Stage-One transfer followed by laboratory procedures for abutment preparation and provisionalization




Common procedures for abutment placement, margin finessing, impressing and placement of restoration






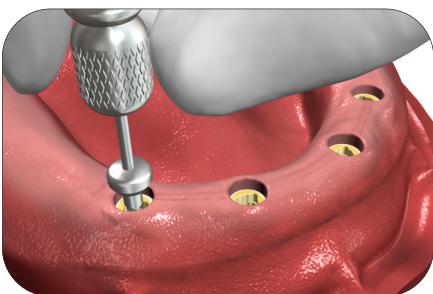
Clinical Procedures for Titanium Temporary Abutment/Fixture-Mount



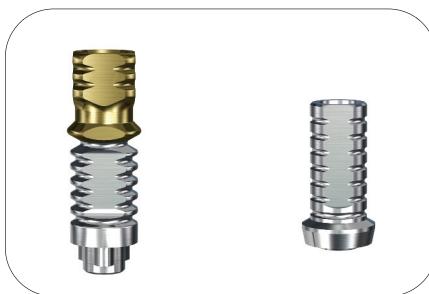
1. Unthread the healing collars from the implant using appropriate driver. The top of the implants are exposed.



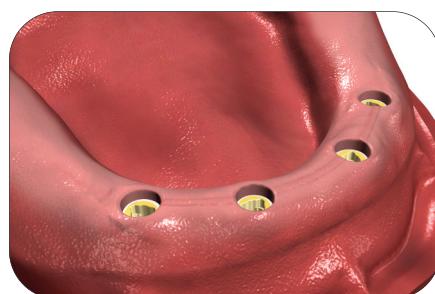2. Use the prepared implant fixture mount as a temporary abutment (reshaping of tri-lobe at apex is essential) or select the Titanium Temporary Abutment matching the implant prosthetic platform.



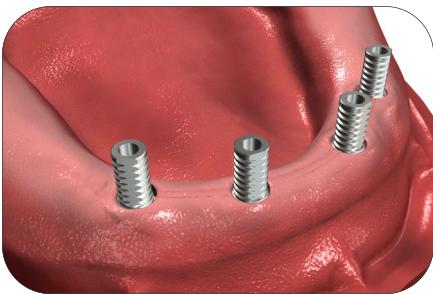
3. Place the fixation screw through the abutment and tighten in place using the appropriate driver. The torque value of the fixation screw is determined by clinician. Standard torque value is 30Ncm.




4. Confirm there is less than 40 degrees of relative divergence between abutments to be splinted together.



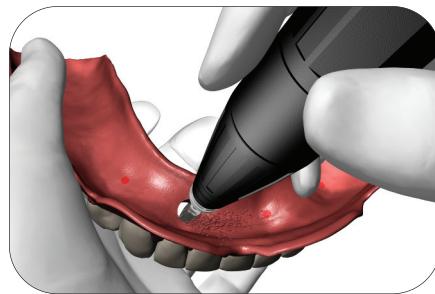5. Mark position of Titanium Temporary Abutments on under surface of denture with indelible marker.



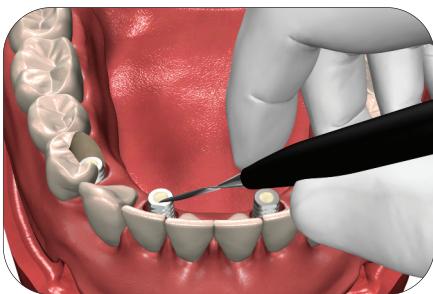6. Use acrylic trimmer to drill holes through denture corresponding to marked positions.



7. Seat denture over Titanium Temporary Abutment to confirm passive placement. Check occlusion with opposing arch. Occlude screw access hole with cotton pellet.



8. Use autopolymerizing acrylic or similar light cure material to attach prosthesis to Titanium Temporary Abutments. Allow to set prior to removal from mouth.



9. Remove prosthesis from the mouth, flip over to expose tissue side. Use additional autopolymerizing resin to fill in any voids around the processed abutments.



10. Optional step: Attach implant analogs to base of denture. Pour a soft tissue working cast to prevent engagement with undercuts. Articulate with opposing arch.




11. Remove denture and trim the prosthesis to allow access for patient cleaning. The prosthesis is now a fixed restoration, removable by clinician only.




12. Attach prosthesis to implants and tighten in place using the appropriate tooling. Occlude screw access holes with material of choice.



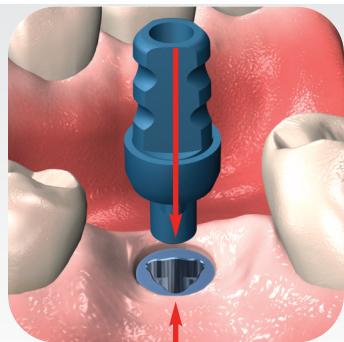
Contoured Abutments with 0, 30, 60 or 90 Degree Tri-Lobe Offset



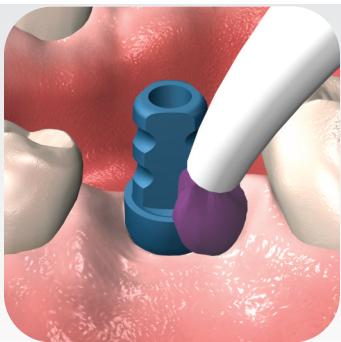
Implant with tri-lobe aligned to labial. Rotate transfer and seat mating tri-lobes.



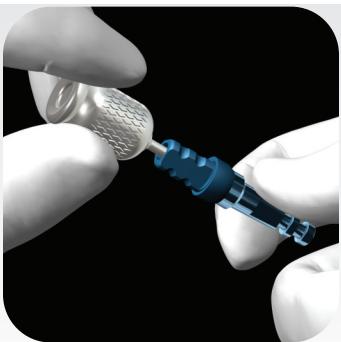Make transfer impression using vinyl polysiloxane material.




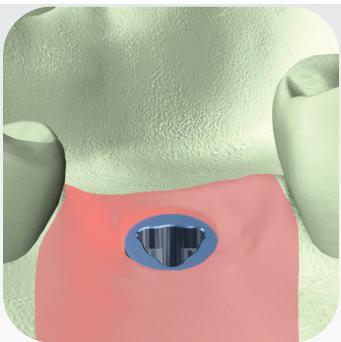Attach transfer to matching colored implant analog. Insert transfer and pour impression.




Evaluate rotational position of tri-lobe for selection of straight and angled contoured abutments

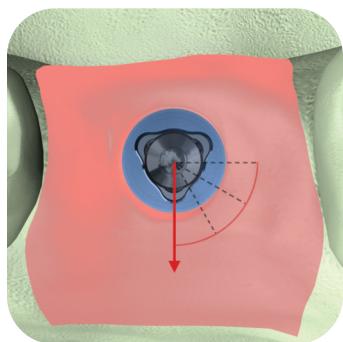0 degrees offset




1. Tri-Lobe aligned perpendicular to labial. (0 degrees offset)

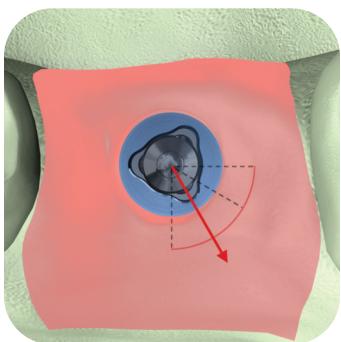30 degrees offset




1. Tri-Lobe aligned 30 degrees counter-clockwise to line perpendicular to labial.

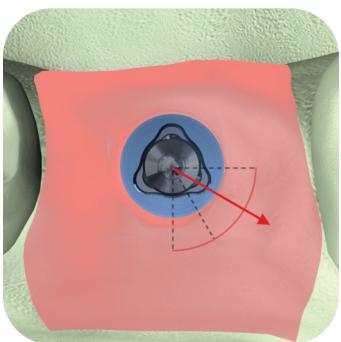60 degrees offset




1. Tri-Lobe aligned 60 degrees counter-clockwise to line perpendicular to labial.

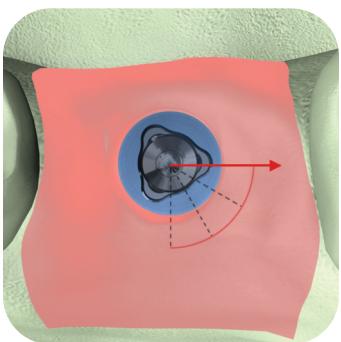
90 degrees offset



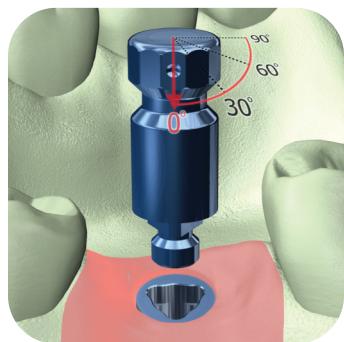
1. Tri-Lobe aligned 90 degrees counter-clockwise to line perpendicular to labial.



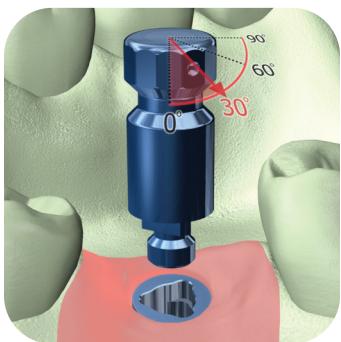
2. Dimple of Alignment Tool is perpendicular to labial: 0 degree offset of implant tri-lobe.



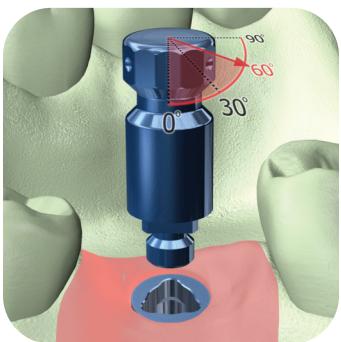
2. Dimple on the tool is 1 facet (30 degrees) counter-clockwise to line perpendicular to labial.



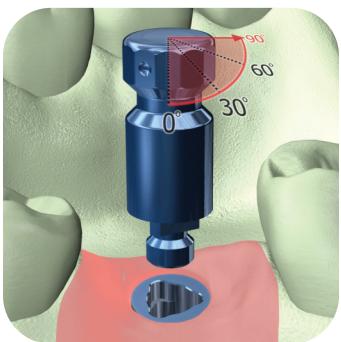
2. Dimple on the tool is 2 facets (60 degrees) counter-clockwise to line perpendicular to labial.



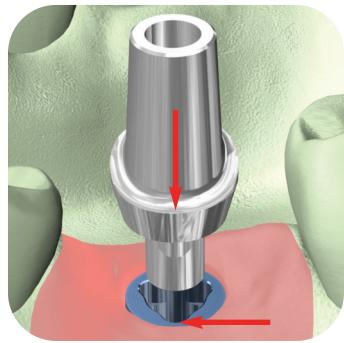
2. Dimple on the tool is 3 facets (90 degrees) counter-clockwise to line perpendicular to labial.



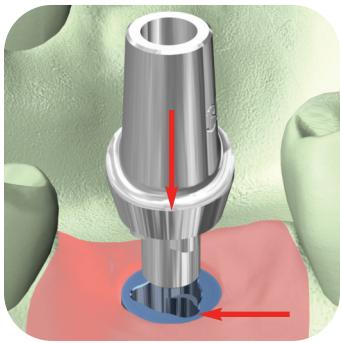
3. Use 0 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



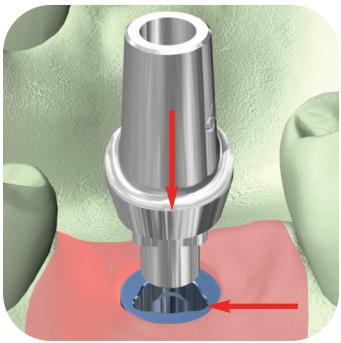3. Use 30 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



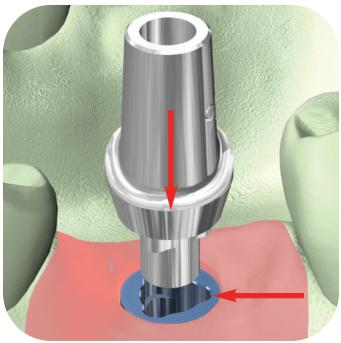3. Use 60 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



3. Use 90 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



Tri-lobe Selection Chart



**Cover Screws, Healing Collars,
Implant Transfers & Analogs**

**Straight Abutments, Comfort
Caps, Transfers & Analogs**

**Straight & Angled
Contoured Abutments**

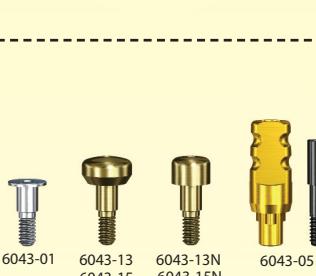
**Lab Abutments &
Engaging & Non-Engaging
Abutments**

**Ball, GPS™ &
Screw-Receiving
Overdenture Abutments**

3.5mmD Platform




4.3mmD Platform




*non-engaging

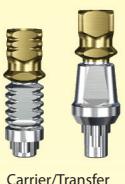5.0mmD Platform




Carrier/Transfer

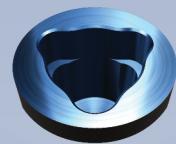

Carrier/Transfer

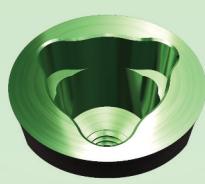6.0mmD Platform



Carrier/Transfer



6043-08 6043-06



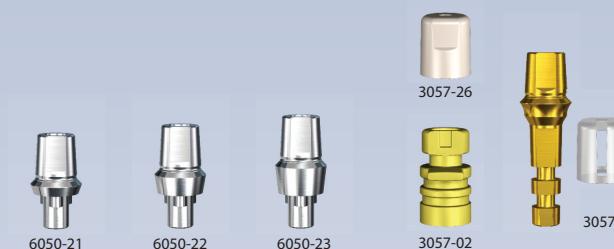
6050-01 6050-13
6050-15 6050-13N
6050-15N 6050-05
6050-08 6050-06



6060-01 6060-13
6060-15 6060-13N
6060-15N 6060-05
6060-08 6060-06



3047-26
3047-02 3047



6050-21 6050-22
6050-23 3057-26
3057-02 3057



6060-21 6060-22
6060-23 3065-26
3065-02 3065



6043-51-15C 6043-51-08C 30/60/90°
6043-52-15C 6043-52-08C 6043-52-30
6043-52-60 6043-52-90



6050-31 6050-31C 30/60/90°
6050-32 6050-32C 6050-32-30
6050-32-60 6050-32-90 6050-51
6050-51-15C 6050-51-08C 30/60/90°
6050-52-30 6050-52-60 6050-52-90



6060-31 6060-32 30/60/90°
6060-32-30 6060-32-60 6060-32-90 6060-51
6060-52 6060-52



6043-47 6043-49



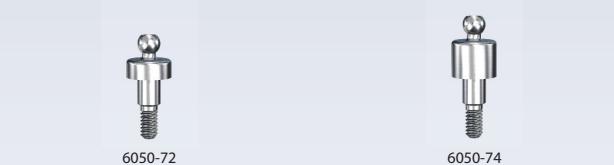
6050-30L 6050-40
6050-41 6050-46 6050-47 6050-49



6060-30L 6060-46 6060-47 6060-49



6043-74



6050-72 6050-74



6060-72 6060-74



6043-64-15 6043-64-30



6050-61 6050-62 6050-64 6050-64-15 6050-64-30



6060-61 6060-62 6060-64



6043-95 6043-94-15 6043-94-30



6050-91 6050-92 6050-93 6050-94 6050-95 6050-94-15 6050-94-30





Healing Collars, Transfers, Analogs & Alignment Tools

Healing Collars

Components include Dual-Grip™ feature for use with 1.25mmD Hex Tool or Nobel Unigrip ScrewDriver


Code	Width	Length	Platform	Price
6035-13	4.6mmD	3mmL	3.5mmD	\$35.00
6035-15	4.6mmD	5mmL	3.5mmD	\$35.00
6043-13	5.6mmD	3mmL	4.3mmD	\$35.00
6043-15	5.6mmD	5mmL	4.3mmD	\$35.00
6050-13	6.1mmD	3mmL	5.0mmD	\$35.00
6050-15	6.1mmD	5mmL	5.0mmD	\$35.00
6060-13	6.6mmD	3mmL	6.0mmD	\$35.00
6060-15	6.6mmD	5mmL	6.0mmD	\$35.00



Narrow Healing Collars

Components include Dual-Grip™ feature for use with 1.25mmD Hex Tool or Nobel Unigrip ScrewDriver

Code	Width	Length	Platform	Price
6035-13N	3.5mmD	3mmL	3.5mmD	\$35.00
6035-15N	3.5mmD	5mmL	3.5mmD	\$35.00
6043-13N	4.3mmD	3mmL	4.3mmD	\$35.00
6043-15N	4.3mmD	5mmL	4.3mmD	\$35.00
6050-13N	5.0mmD	3mmL	5.0mmD	\$35.00
6050-15N	5.0mmD	5mmL	5.0mmD	\$35.00
6060-13N	6.0mmD	3mmL	6.0mmD	\$35.00
6060-15N	6.0mmD	5mmL	6.0mmD	\$35.00

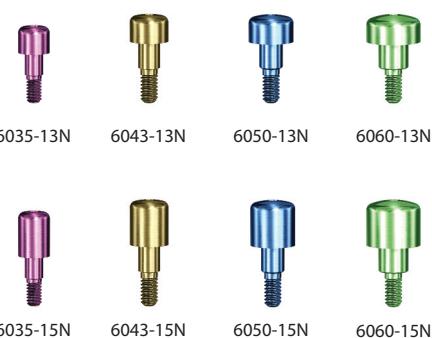

Closed-Tray Implant Transfers

Includes Dual-Grip™ Fixation Screw

Code	Platform	Replacement Screw - 2ea	Price
6035-05	3.5mmD	1000-28	\$35.00
6043-05	4.3mmD	1000-27	\$35.00
6050-05	5.0mmD	1000-27	\$35.00
6060-05	6.0mmD	1000-27	\$35.00



Open-Tray Implant Transfers


Includes Long Fixation Screw
[Requires 1.25mmD Hex Tool]

Code	Platform	Replacement Screw - 2ea	Price
6035-08	3.5mmD	1000-12	\$35.00
6043-08	4.3mmD	1000-11	\$35.00
6050-08	5.0mmD	1000-11	\$35.00
6060-08	6.0mmD	1000-11	\$35.00



Implant Analogs

Code	Platform	Price
6035-06	3.5mmD	\$20.00
6043-06	4.3mmD	\$20.00
6050-06	5.0mmD	\$20.00
6060-06	6.0mmD	\$20.00

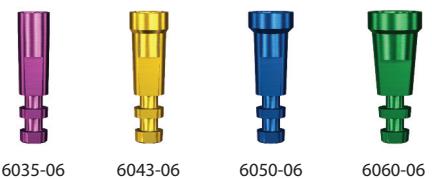

Lab Abutment Alignment Tools

Code	Platform	Price
6035-07	3.5mmD	\$40.00
6043-07	4.3mmD	\$40.00
6050-07	5.0mmD	\$40.00
6060-07	6.0mmD	\$40.00






Straight Abutments, Snap-on Transfers & Lab Analogs

Straight Abutments

Includes Transfer, Comfort Cap & Fixation Screw with Dual-Grip™

Code	Platform	Collar Heights	Use Analog	Replacement Screw	Price
6035-21	3.5mmD	1mm	3037	1000-23	\$100.00
6035-22	3.5mmD	2mm	3037	1000-23	\$100.00
6035-23	3.5mmD	3mm	3037	1000-23	\$100.00
6043-21	4.3mmD	1mm	3047	1000-22	\$100.00
6043-22	4.3mmD	2mm	3047	1000-22	\$100.00
6043-23	4.3mmD	3mm	3047	1000-22	\$100.00
6050-21	5.0mmD	1mm	3057	1000-22	\$100.00
6050-22	5.0mmD	2mm	3057	1000-22	\$100.00
6050-23	5.0mmD	3mm	3057	1000-22	\$100.00
6060-21	6.0mmD	1mm	3065	1000-22	\$100.00
6060-22	6.0mmD	2mm	3065	1000-22	\$100.00
6060-23	6.0mmD	3mm	3065	1000-22	\$100.00

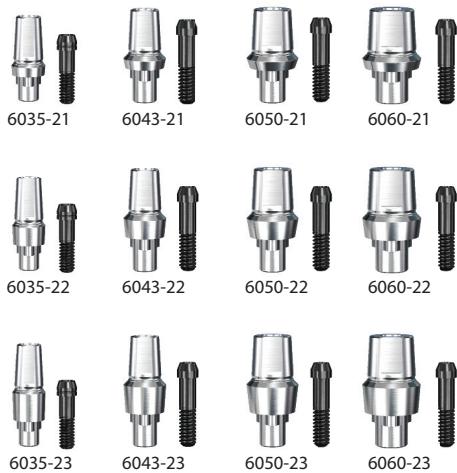

Straight Abutment Transfers

Code	Abutment Width	Price (includes two pieces)
3037-02	3.7mmD	\$30.00
3047-02	4.7mmD	\$30.00
3057-02	5.7mmD	\$30.00
3065-02	6.5mmD	\$30.00



Straight Abutment Comfort Caps

Code	Abutment Width	Price (includes two pieces)
3037-26	3.7mmD	\$30.00
3047-26	4.7mmD	\$30.00
3057-26	5.7mmD	\$30.00
3065-26	6.5mmD	\$30.00



Straight Abutment Analogs

Includes Castable Coping

Code	Abutment Width	Price
3037	3.7mmD	\$35.00
3047	4.7mmD	\$35.00
3057	5.7mmD	\$35.00
3065	6.5mmD	\$35.00



Straight Abutment Castable Copings

Code	Abutment Width	Price (includes two pieces)
3037-28	3.7mmD	\$30.00
3047-28	4.7mmD	\$30.00
3057-28	5.7mmD	\$30.00
3065-28	6.5mmD	\$30.00





Straight & Angled Contoured, Titanium & Zirconia Abutments

Straight Titanium Contoured Abutments

Includes Dual-Grip™ Fixation Screw

Code	Platform	Collar Heights	Replacement Screw	Price
6035-31	3.5mmD	1mm	1000-23	\$100.00
6035-32	3.5mmD	2mm	1000-23	\$100.00
6043-31	4.3mmD	1mm	1000-22	\$100.00
6043-32	4.3mmD	2mm	1000-22	\$100.00
6050-31	5.0mmD	1mm	1000-22	\$100.00
6050-32	5.0mmD	2mm	1000-22	\$100.00
6060-31	6.0mmD	1mm	1000-22	\$100.00
6060-32	6.0mmD	2mm	1000-22	\$100.00



15° Angled Titanium Contoured Abutments

Includes Dual-Grip™ Fixation Screw

Code	Platform	Collar Heights	Replacement Screw	Price
6035-51	3.5mmD	.75mm	1000-23	\$100.00
6035-52	3.5mmD	2mm	1000-23	\$100.00
6043-51	4.3mmD	.75mm	1000-22	\$100.00
6043-52	4.3mmD	2mm	1000-22	\$100.00
6050-51	5.0mmD	.75mm	1000-22	\$100.00
6050-52	5.0mmD	2mm	1000-22	\$100.00
6060-51	6.0mmD	.75mm	1000-22	\$100.00
6060-52	6.0mmD	2mm	1000-22	\$100.00



0°, 8° & 15° Zirconia Contoured Abutments

Includes Dual-Grip™ Fixation Screw

Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
6035-31C	3.5mmD	0°	1mm	1000-23	\$120.00
6035-32C	3.5mmD	0°	2mm	1000-23	\$120.00
6035-51-08C	3.5mmD	8°	1mm	1000-23	\$120.00
6035-52-08C	3.5mmD	8°	2mm	1000-23	\$120.00
6035-51-15C	3.5mmD	15°	1mm	1000-23	\$120.00
6035-52-15C	3.5mmD	15°	2mm	1000-23	\$120.00
6035-30CB	3.5mmD	-	.5mm	1000-23	\$60.00



Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
6043-31C	4.3mmD	0°	1mm	1000-22	\$120.00
6043-32C	4.3mmD	0°	2mm	1000-22	\$120.00
6043-51-08C	4.3mmD	8°	1mm	1000-22	\$120.00
6043-52-08C	4.3mmD	8°	2mm	1000-22	\$120.00
6043-51-15C	4.3mmD	15°	1mm	1000-22	\$120.00
6043-52-15C	4.3mmD	15°	2mm	1000-22	\$120.00
6043-30CB	4.3mmD	-	.5mm	1000-22	\$60.00



Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
6050-31C	5.0mmD	0°	1mm	1000-22	\$120.00
6050-32C	5.0mmD	0°	2mm	1000-22	\$120.00
6050-51-08C	5.0mmD	8°	1mm	1000-22	\$120.00
6050-52-08C	5.0mmD	8°	2mm	1000-22	\$120.00
6050-51-15C	5.0mmD	15°	1mm	1000-22	\$120.00
6050-52-15C	5.0mmD	15°	2mm	1000-22	\$120.00
6050-30CB	5.0mmD	-	.5mm	1000-22	\$60.00





Angled, Contoured & Lab Titanium Abutments

Straight Contoured Abutments with 30, 60 & 90 Degree offset Tri-Lobe

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-32-30	2mm	3.5mmD	1000-23	\$100.00
6035-32-60	2mm	3.5mmD	1000-23	\$100.00
6035-32-90	2mm	3.5mmD	1000-23	\$100.00
6043-32-30	2mm	4.3mmD	1000-22	\$100.00
6043-32-60	2mm	4.3mmD	1000-22	\$100.00
6043-32-90	2mm	4.3mmD	1000-22	\$100.00
6050-32-30	2mm	5.0mmD	1000-22	\$100.00
6050-32-60	2mm	5.0mmD	1000-22	\$100.00
6050-32-90	2mm	5.0mmD	1000-22	\$100.00
6060-32-30	2mm	6.0mmD	1000-22	\$100.00
6060-32-60	2mm	6.0mmD	1000-22	\$100.00
6060-32-90	2mm	6.0mmD	1000-22	\$100.00



Angled Contoured Abutments with 30, 60 & 90 Degree offset Tri-Lobe

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-52-30	2mm	3.5mmD	1000-23	\$100.00
6035-52-60	2mm	3.5mmD	1000-23	\$100.00
6035-52-90	2mm	3.5mmD	1000-23	\$100.00
6043-52-30	2mm	4.3mmD	1000-22	\$100.00
6043-52-60	2mm	4.3mmD	1000-22	\$100.00
6043-52-90	2mm	4.3mmD	1000-22	\$100.00
6050-52-30	2mm	5.0mmD	1000-22	\$100.00
6050-52-60	2mm	5.0mmD	1000-22	\$100.00
6050-52-90	2mm	5.0mmD	1000-22	\$100.00



Straight Full Contour Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-30L	-	3.5mmD	1000-23	\$85.00
6043-30L	-	4.3mmD	1000-22	\$85.00
6050-30L	-	5.0mmD	1000-22	\$85.00
6060-30L	-	6.0mmD	1000-22	\$85.00



Replacement Fixation Screws

*Components include Dual-Grip™ feature for use with
1.25mmD Hex Tool or Nobel Unigrip ScrewDriver*

Code	Platform	Price (includes two pieces)
1000-23	3.5mmD only (short)	\$30.00
1000-28	3.5mmD only (long)	\$30.00
1000-22	4.3, 5.0 and 6.0mmD (short)	\$30.00
1000-27	4.3, 5.0 and 6.0mmD (long)	\$30.00





Engaging & Non-Engaging Abutments

Gold Engaging Cast-to Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-40	1mm	3.5mmD	1000-23	\$120.00
6043-40	1mm	4.3mmD	1000-22	\$120.00
6050-40	1mm	5.0mmD	1000-22	\$120.00



Gold Non-Engaging Cast-to Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-41	1mm	3.5mmD	1000-23	\$120.00
6043-41	1mm	4.3mmD	1000-22	\$120.00
6050-41	1mm	5.0mmD	1000-22	\$120.00



Plastic Non-Engaging Castable Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-46	-	3.5mmD	1000-23	\$35.00
6043-46	-	4.3mmD	1000-22	\$35.00
6050-46	-	5.0mmD	1000-22	\$35.00
6060-46	-	6.0mmD	1000-22	\$35.00



Plastic Engaging Temporary Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-47	-	3.5mmD	1000-28	\$40.00
6043-47	-	4.3mmD	1000-27	\$40.00
6050-47	-	5.0mmD	1000-27	\$40.00
6060-47	-	6.0mmD	1000-27	\$40.00



Titanium Non-Engaging Temporary Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Replacement Screw	Price
6035-49	-	3.5mmD	1000-23	\$40.00
6043-49	-	4.3mmD	1000-22	\$40.00
6050-49	-	5.0mmD	1000-22	\$40.00
6060-49	-	6.0mmD	1000-22	\$40.00





Selection Chart: Overdenture & Screw-Retained Restorations

3.5mmD Platform



4.3mmD Platform



5.0mmD Platform



6.0mmD Platform



Overdenture Abutments (*Includes Transfer and Comfort Cap)

Overdenture Abutments



Overdenture Abutment Transfers, Analog & Comfort Cap

OD Abutment Comfort Cap Transfers and Analog

*4000-13

*1000-68



Overdenture Copings

Overdenture Abutment Components



1000-64



1000-65



1000-66



1000-67



1000-69

Ball Abutments (**Includes Transfer and Comfort Cap)

Ball Abutments



Ball Attachments

Ball Abutment Components



**1000-87



**1000-84




1000-85



1000-81



1000-82



1000-88



1000-83



1000-86

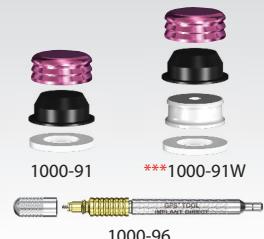
GPS™ Abutments (**Includes Transfer, Abutment Housing Kit and Comfort Cap)

GPS™ Abutments



GPS™ Attachments

GPS™ Abutment Components






Screw-Receiving & Ball Abutments

Overdenture Abutments

*Includes Transfer & Comfort Cap
[Requires 1.25mmD Hex Tool]*

Code	Collar Heights	Platform	Price
6035-61	1mm	3.5mmD	\$100.00
6035-62	2mm	3.5mmD	\$100.00
6035-64	4mm	3.5mmD	\$100.00
6043-61	1mm	4.3mmD	\$100.00
6043-62	2mm	4.3mmD	\$100.00
6043-64	4mm	4.3mmD	\$100.00
6050-61	1mm	5.0mmD	\$100.00
6050-62	2mm	5.0mmD	\$100.00
6050-64	4mm	5.0mmD	\$100.00
6060-61	1mm	6.0mmD	\$100.00
6060-62	2mm	6.0mmD	\$100.00
6060-64	4mm	6.0mmD	\$100.00

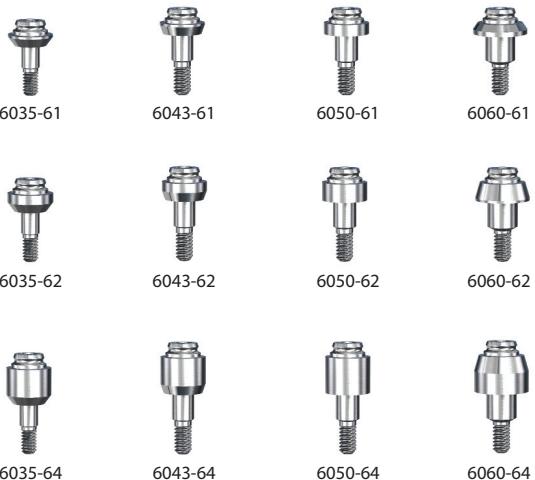

Angled Overdenture Abutments

Includes Dual-Grip™ Abutment Screw, Transfer & Comfort Cap [Screw Top requires 2.5mmD Hex Tool]


Code	Description	Platform	Price
6035-64-15	15° Angled Abutment	3.5mmD	\$120.00
6043-64-15	15° Angled Abutment	4.3mmD	\$120.00
6050-64-15	15° Angled Abutment	5.0mmD	\$120.00
6035-64-30	30° Angled Abutment	3.5mmD	\$120.00
6043-64-30	30° Angled Abutment	4.3mmD	\$120.00
6050-64-30	30° Angled Abutment	5.0mmD	\$120.00



Ball Abutments

*Includes Transfer & Comfort Cap
[Requires 1.25mmD Hex Tool]*

Code	Collar Heights	Platform	Price
6035-72	2mm	3.5mmD	\$90.00
6035-74	4mm	3.5mmD	\$90.00
6043-72	2mm	4.3mmD	\$90.00
6043-74	4mm	4.3mmD	\$90.00
6050-72	2mm	5.0mmD	\$90.00
6050-74	4mm	5.0mmD	\$90.00
6060-72	2mm	6.0mmD	\$90.00
6060-74	4mm	6.0mmD	\$90.00

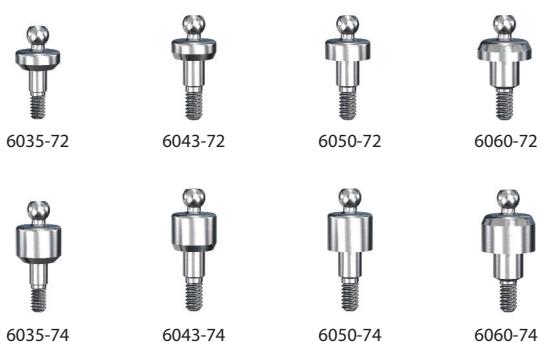



GPS™ Abutments & Angled Base Components

Angled GPS Abutments

Includes Abutment Screw, Transfer, Comfort Cap, GPS Abutment Housing Kit & White Male Liner [GPS Top requires 1.25mmD Hex Tool]

Code	Description	Platform	Price
6035-94-15	15° Angled Abutment	3.5mmD	\$120.00
6043-94-15	15° Angled Abutment	4.3mmD	\$120.00
6050-94-15	15° Angled Abutment	5.0mmD	\$120.00
6035-94-30	30° Angled Abutment	3.5mmD	\$120.00
6043-94-30	30° Angled Abutment	4.3mmD	\$120.00
6050-94-30	30° Angled Abutment	5.0mmD	\$120.00



GPS Abutments

Includes Transfer, Comfort Cap, GPS Abutment Housing Kit & White Male Liner [use with 1.25mmD Hex Tool or Nobel Unigrip ScrewDriver]

Code	Collar Heights	Platform	Price
6035-91	1mm	3.5mmD	\$100.00
6035-92	2mm	3.5mmD	\$100.00
6035-93	3mm	3.5mmD	\$100.00
6035-94	4mm	3.5mmD	\$100.00
6035-95	5mm	3.5mmD	\$100.00
6043-91	1mm	4.3mmD	\$100.00
6043-92	2mm	4.3mmD	\$100.00
6043-93	3mm	4.3mmD	\$100.00
6043-94	4mm	4.3mmD	\$100.00
6043-95	5mm	4.3mmD	\$100.00
6050-91	1mm	5.0mmD	\$100.00
6050-92	2mm	5.0mmD	\$100.00
6050-93	3mm	5.0mmD	\$100.00
6050-94	4mm	5.0mmD	\$100.00
6050-95	5mm	5.0mmD	\$100.00



Replacement Tops for Angled Abutments

Includes Transfer & Comfort Cap

Code	Description	Price
1000-60	Replacement Screw Top for Angled Base [Requires 2.5mmD Hex Tool]	\$35.00
1000-70	Replacement Ball Top for Angled Base [Requires 1.25mmD Hex Tool]	\$35.00
1000-90	Replacement GPS Top for Angled Base [Requires 1.25mmD Hex Tool]	\$35.00





Ancillary Components for Overdenture & Screw-Retained Restorations

Overdenture Abutment Comfort Cap, Transfers & Analog

Screws require 1.25mmD Hex Tool


Code	Platform	Description	Price
1000-61	5.0mmD	Aluminum Analog	\$30.00
1000-62	5.0mmD	Open-Tray Transfer plus Fixation Screw	\$35.00
1000-63	5.0mmD	Closed-Tray Transfer plus Fixation Screw	\$35.00
1000-68	5.0mmD	Plastic Abutment Transfer - 2ea	\$30.00
4000-13	5.0mmD	Plastic Comfort Cap - 2ea	\$30.00
1000-14	-	Replacement Open-Tray Fixation Screw - 2ea	\$30.00
1000-15	-	Replacement Closed-Tray Fixation Screw - 2ea	\$30.00



Overdenture Abutment Components

Includes Fixation Screw [1000-05], optional Dual-Grip™ Fixation Screw [1000-25]

Code	Platform	Description	Price
1000-64	5.0mmD	Gold Coping	\$80.00
1000-65	5.0mmD	Titanium Coping	\$40.00
1000-66	5.0mmD	Plastic Coping	\$35.00
1000-67	5.0mmD	Gold & Plastic Coping	\$80.00
1000-69	5.0mmD	Titanium Temporary Coping	\$35.00
1000-05	-	Coping Fixation Screw - 2ea	\$30.00
1000-25	-	Dual-Grip™ Coping Fixation Screw - 2ea	\$30.00

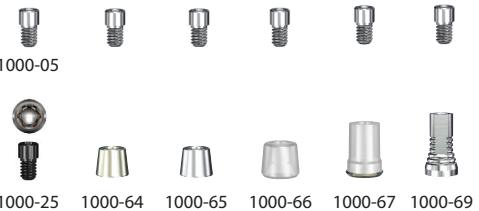

Ball Abutment Components

Code	Description	Price
1000-81	Ball Abutment Housing - 2ea	\$30.00
1000-82	Nylon Liner - 2ea	\$15.00
1000-83	Nylon Liner Transfer - 2ea	\$15.00
1000-84	Ball Abutment Transfer - 2ea	\$30.00
1000-85	Ball Abutment Analog	\$20.00
1000-86	Ball Abutment Tools	\$65.00
1000-87	Ball Abutment Comfort Cap - 2ea	\$30.00
1000-88	Ball Abutment Housing & Nylon Liner	\$22.00



GPS Abutment Transfer, Analogs & Alignment Tools

Code	Description	Price
1000-94	Impression Coping - 2ea	\$30.00
1000-95	Abutment Analog 4mmD - 2ea	\$40.00
1000-97	Abutment Parallel Post - 4ea	\$16.00
4500-26	Comfort Cap - 2ea	\$10.00
1000-99	Plastic Abutment Angle Measure Guide	\$5.00



GPS Abutment Components

Code	Retention Force	Description	Price
1000-92W	5.0lbs	Internal Replacement Male White - 4ea	\$22.00
1000-92P	3.0lbs	Internal Replacement Male Pink - 4ea	\$22.00
1000-92B	1.5lbs	Internal Replacement Male Blue - 4ea	\$22.00
1000-92C	4.5lbs	External Replacement Male Clear - 4ea	\$16.00
1000-92G	3.0lbs	External Replacement Male Green - 4ea	\$16.00
1000-92E	2.0lbs	External Replacement Male Orange - 4ea	\$16.00
1000-92R	1.0lb	External Replacement Male Red - 4ea	\$16.00
1000-91	-	GPS Abutment Housing Kit - 1ea	\$10.00
1000-91W	-	GPS Housing Kit w/ White Internal Male- 1ea	\$15.00
1000-93	-	Black Processing Male - 4ea	\$16.00
1000-98SI	Assorted	Internal Processing Cap/Male Kit - 1ea	\$26.00
1000-98S	Assorted	External Processing Cap/Male Kit - 1ea	\$20.00



GPS Instrumentation

Code	Description	Price
1000-96	Core Tool - Male Removal & Seating Tool	\$75.00





Surgical & Prosthetic Instruments


Tapered Hex Tools

Code	Description	Price
THT1.25	1.25mmD Tapered Hex Tool-Long	\$40.00
THT1.25S	1.25mmD Tapered Hex Tool-Short	\$40.00



Straight Hex Tools

Code	Description	Price
HT1.25L	1.25mmD Straight Hex Tool Extra Long (For use with Torque Wrench)	\$40.00
HT1.25	1.25mmD Straight Hex Tool Long (For use with Torque Wrench)	\$40.00
HT1.25S	1.25mmD Straight Hex Tool Short (For use with Torque Wrench)	\$40.00
HT2.0	2.0mmD Straight Hex Tool Long	\$40.00
HT2.0S	2.0mmD Straight Hex Tool Short	\$40.00
HT2.5	2.5mmD Straight Hex Tool-Long (with retentive O-ring)	\$40.00
HT2.5S	2.5mmD Straight Hex Tool-Short (with retentive O-ring)	\$40.00
HT3.0	3.0mmD Straight Hex Tool-Long	\$40.00
HT3.0S	3.0mmD Straight Hex Tool-Short	\$40.00

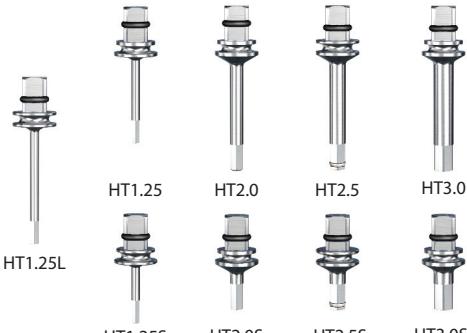



Hex Drills

Code	Description	Price
HD1.25	1.25mmD Hex Drill	\$40.00
HD2.0	2.0mmD Hex Drill	\$40.00
HD2.5	2.5mmD Hex Drill (with retentive O-ring)	\$40.00
HD3.0	3.0mmD Hex Drill	\$40.00




Ancillary Instruments

Code	Description	Price
R	Ratchet	\$125.00
TW30	30Ncm Torque Wrench	\$125.00
TW35	35Ncm Torque Wrench	\$125.00
SSH2.5	2.5mmD Stainless Steel Handle (with retentive O-ring)	\$100.00
FMD	Fixture Mount Drill (with retentive O-ring)	\$80.00
PAR	Paralleling Tool 2.8/2.3	\$15.00
LDRILL	Locator Drill (2.3mmD)	\$30.00
DE	Drill Extender	\$40.00
TWM	Multi-Setting Torque Wrench (15, 20, 25, 30, 35, 40, 45, 50 & 60Ncm)	\$395.00
ABH	Analog Holder	\$85.00

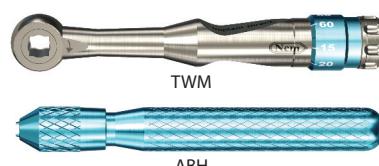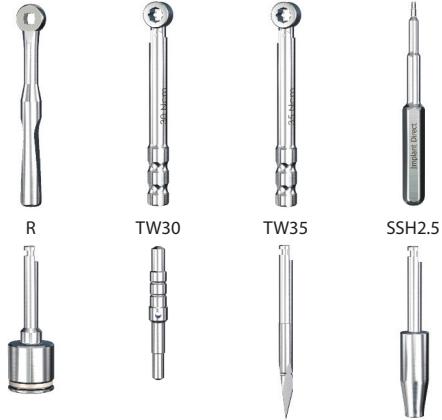

Surgical Trays

Code	Description	Price
ST	Empty Surgical Kit for RePlus & ReActive System Implants	\$250.00
BST	Basic Surgical Kit for RePlus & ReActive System Implants	\$1265.00
CST	Surgical Kit for RePlus & ReActive System Implants	\$2565.00
RST	Empty Surgical Kit for RePlant System Implants	\$250.00
RBST	Basic Surgical Kit for RePlant System Implants	\$1435.00
RCST	Surgical Kit for RePlant System Implants	\$2325.00

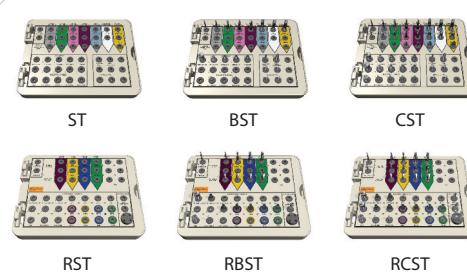



Surgical Drills: RePlus & ReActive Implants

Pilot Drills

Code	Description	Price
SD2.3	2.3/2.0mmD Step Drill with 6-16 markings	\$60.00
SD2.3S	2.3/2.0mmD Step Drill with 6-10 markings	\$60.00
SD2.8	2.8/2.3mmD Step Drill with 6-16 markings	\$60.00
SD2.8S	2.8/2.3mmD Step Drill with 6-10 markings	\$60.00

SD2.3



SD2.3S



SD2.8



SD2.8S



3.7mmD Implants

Code	Description	Price
SD3.4	3.4/2.8mmD Step Drill with 6-16 markings	\$60.00
SD3.4S	3.4/2.8mmD Step Drill with 6-10 markings	\$60.00
CBD3.7	Cortical Bone Drill for 3.7mmD Implants	\$60.00

SD3.4



SD3.4S



CBD3.7



4.2mmD Implants

Code	Description	Price
SD3.8	3.8/3.4mmD Step Drill with 6-16 markings	\$60.00
SD3.8S	3.8/3.4mmD Step Drill with 6-10 markings	\$60.00
CBD4.2	Cortical Bone Drill for 4.2mmD Implants	\$60.00

SD3.8



SD3.8S



CBD4.2



4.7mmD Implants

Code	Description	Price
SD4.4	4.4/3.8mmD Step Drill with 6-16 markings	\$60.00
SD4.4S	4.4/3.8mmD Step Drill with 6-10 markings	\$60.00
CBD4.7	Cortical Bone Drill for 4.7mmD Implants	\$60.00

SD4.4



SD4.4S



CBD4.7



Optional Drills

Code	Description	Price
SD4.8	4.8/4.4mmD Step Drill with 6-16 markings	\$60.00
SD4.8S	4.8/4.4mmD Step Drill with 6-10 markings	\$60.00

SD4.8



SD4.8S



5.7mmD Implants

Code	Description	Price
D5.1	5.1mmD Straight Drill with 6-16 markings	\$60.00
D5.1S	5.1mmD Straight Drill with 6-10 markings	\$60.00
SD5.4	5.4/4.8mmD Step Drill with 6-16 markings	\$60.00
SD5.4S	5.4/4.8mmD Step Drill with 6-10 markings	\$60.00
CBD5.7	Cortical Bone Drill for 5.7mmD Implants	\$60.00

D5.1



D5.1S



SD5.4



SD5.4S



CBD5.7



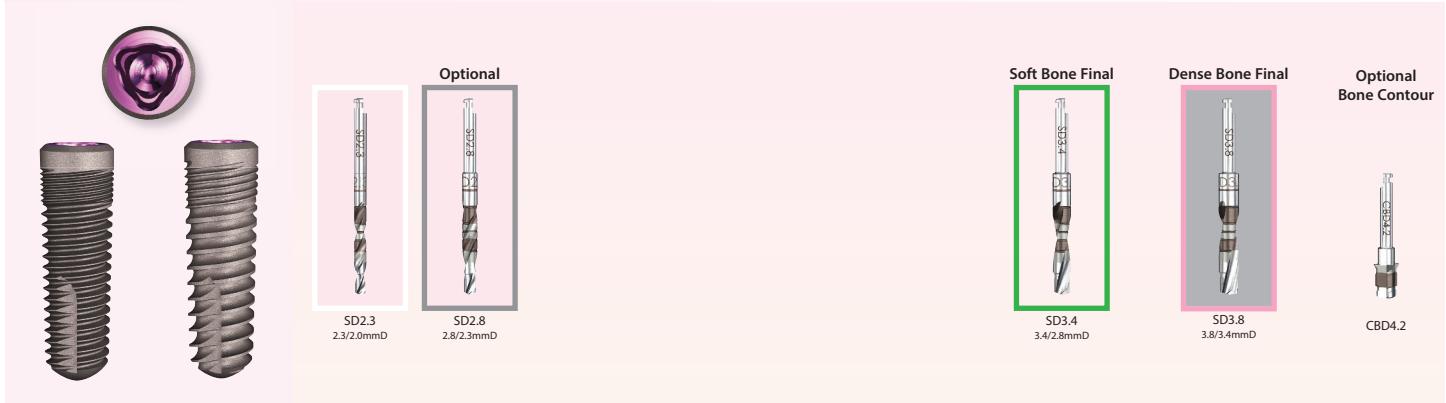


RePlus & ReActive Drill Sequence


3.7mmD RePlus & ReActive Tapered Implant




4.2mmD RePlus & ReActive Tapered Implant



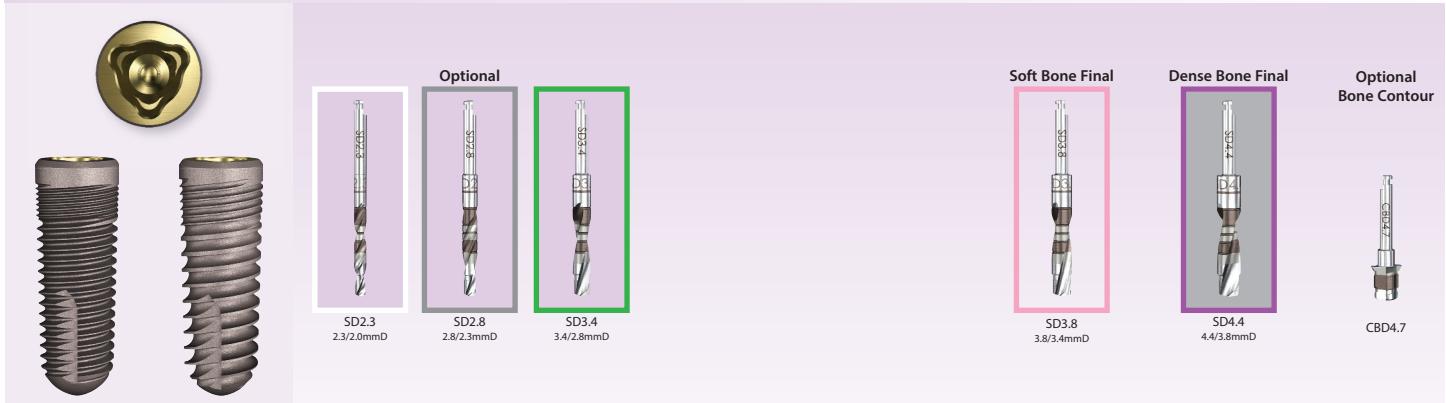4.7mmD RePlus & ReActive Tapered Implant




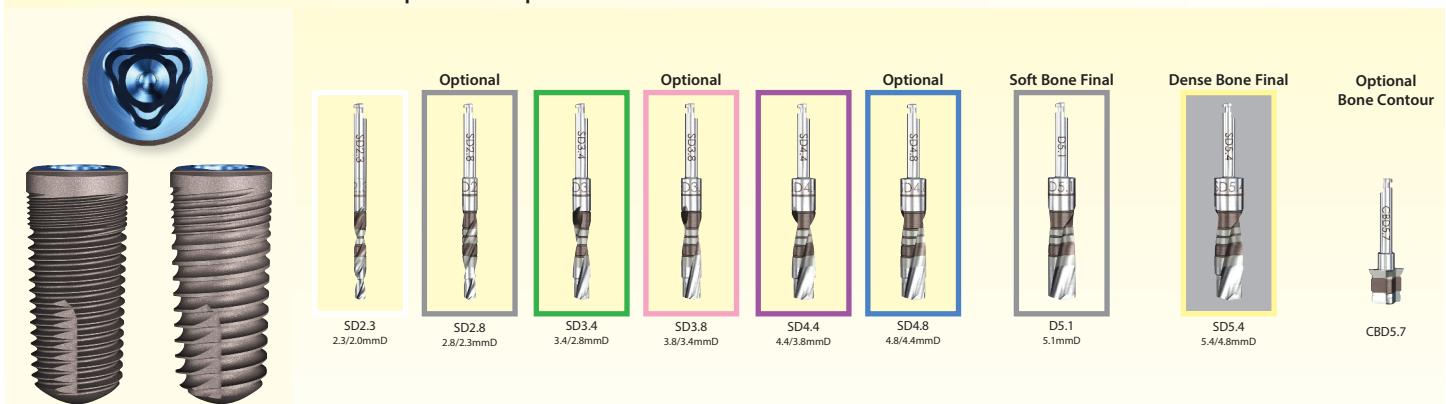5.7mmD RePlus & ReActive Tapered Implant





Complete Surgical Tray Contents

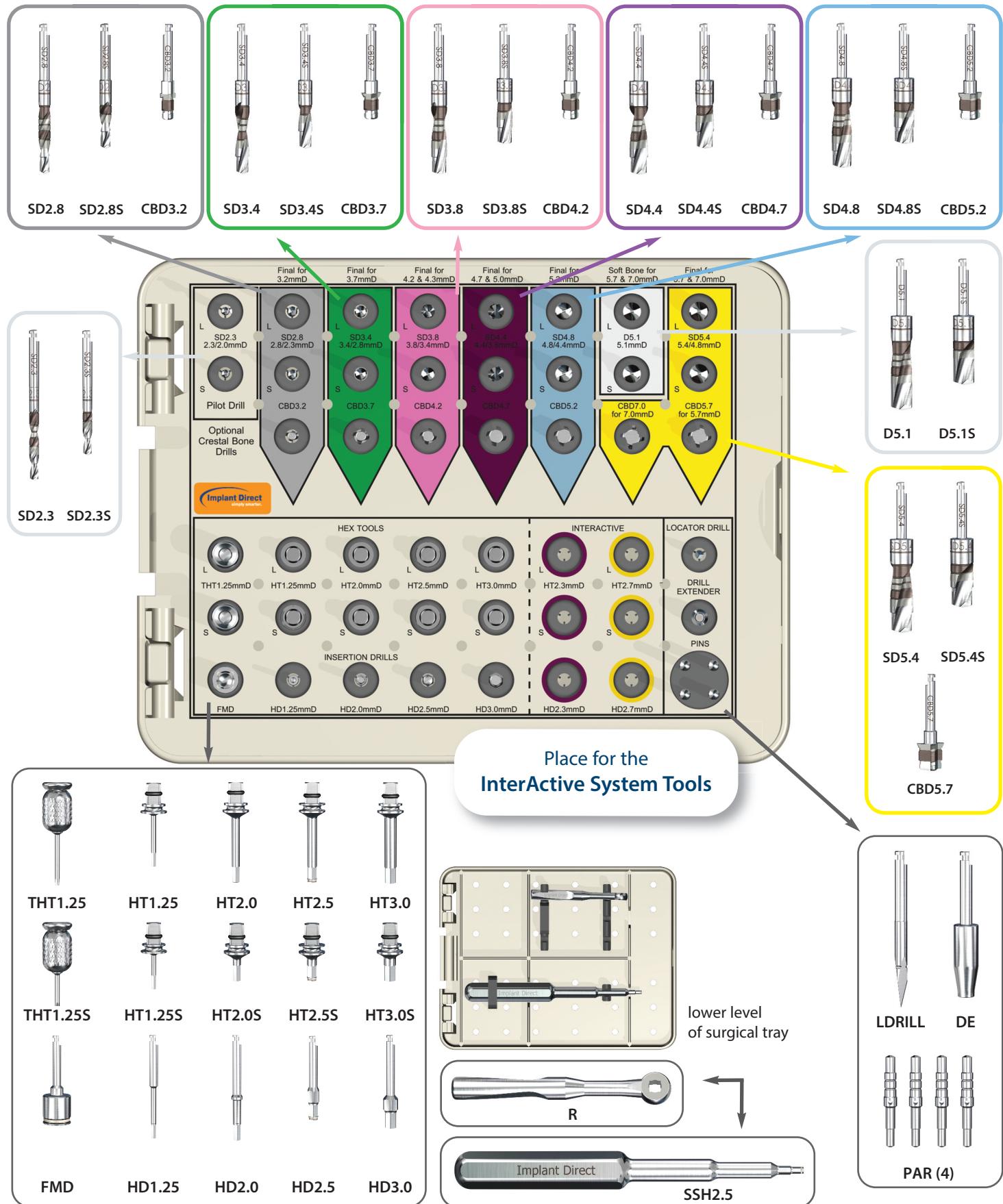





Surgical Tray



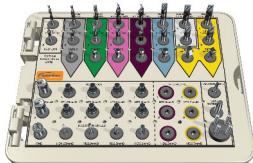Lid



Inside

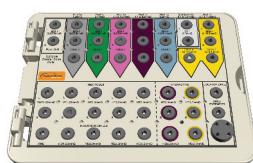Empty
Tray

Surgical Tray Contents			Complete Kit	Basic Kit	Conversion Kit
Code	Description	Price			
ST	Empty Surgical Tray	\$250.00	X	X	X
SD2.3	2.3/2.0mmD Step Drill-Long	\$60.00	X	X	
SD2.3S	2.3/2.0mmD Step Drill-Short	\$60.00	X		
SD2.8	2.8/2.3mmD Step Drill-Long	\$60.00	X	X	
SD2.8S	2.8/2.3mmD Step Drill-Short	\$60.00	X		
CBD3.2	3.2mmD Cortical Drill	\$60.00	X		X
SD3.4	3.4/2.8mmD Step Drill-Long	\$60.00	X	X	
SD3.4S	3.4/2.8mmD Step Drill-Short	\$60.00	X		
CBD3.7	3.7mmD Cortical Drill	\$60.00	X		X
SD3.8	3.8/3.4mmD Step Drill-Long	\$60.00	X	X	X
SD3.8S	3.8/3.4mmD Step Drill-Short	\$60.00	X		X
CBD4.2	4.2mmD Cortical Drill	\$60.00	X		X
SD4.4	4.4/3.8mmD Step Drill-Long	\$60.00	X	X	
SD4.4S	4.4/3.8mmD Step Drill-Short	\$60.00	X		
CBD4.7	4.7mmD Cortical Drill	\$60.00	X		X
SD4.8	4.8/4.4mmD Step Drill-Long	\$60.00	X	X	X
SD4.8S	4.8/4.4mmD Step Drill-Short	\$60.00	X		X
CBD5.2	5.2mmD Cortical Drill	\$60.00	X		X
D5.1	5.1mmD Straight Drill-Long	\$60.00	X	X	
D5.1S	5.1mmD Straight Drill-Short	\$60.00	X		
SD5.4	5.4/4.8mmD Step Drill-Long	\$60.00	X	X	
SD5.4S	5.4/4.8mmD Step Drill-Short	\$60.00	X		
CBD5.7	5.7mmD Cortical Drill	\$60.00	X		X
Hex Tools					
THT1.25	1.25mmD Tapered Hex Tool-Long	\$40.00	X	X	
THT1.25S	1.25mmD Tapered Hex Tool-Short	\$40.00	X		
HT1.25	1.25mmD Straight Hex Tool-Long	\$40.00	X	X	
HT1.25S	1.25mmD Straight Hex Tool-Short	\$40.00	X		
HT2.0	2.0mmD Hex Tool-Long	\$40.00	X	X	
HT2.0S	2.0mmD Hex Tool-Short	\$40.00	X		
HT2.5	2.5mmD Hex Tool-Long	\$40.00	X	X	
HT2.5S	2.5mmD Hex Tool-Short	\$40.00	X		
HT3.0	3.0mmD Hex Tool-Long	\$40.00	X	X	
HT3.0S	3.0mmD Hex Tool-Short	\$40.00	X		
HD1.25	1.25mmD Hex Drill	\$40.00	X	X	
HD2.0	2.0mmD Hex Drill	\$40.00	X		X
HD2.5	2.5mmD Hex Drill	\$40.00	X	X	
HD3.0	3.0mmD Hex Drill	\$40.00	X		X
Ancillary Instrumentation					
FMD	Fixture Mount Drill	\$80.00	X		
LDRILL	Locator Drill	\$30.00	X	X	
DE	Drill Extender	\$40.00	X	X	
R	Ratchet	\$125.00	X	X	
PAR (x4)	Paralleling Tool 2.8/2.3	\$60.00	X	X	
SSH2.5	2.5mmD Stainless Steel Handle	\$100.00	X		



Surgical Drills: RePlant® & NobelReplace™ Tapered Groovy Implants

Pilot Drills for All Implants

Code	Description	Price
SD2.3	2.3/2.0mmD Step Drill with 6-16 markings	\$60.00
SD2.3S	2.3/2.0mmD Step Drill with 6-10 markings	\$60.00

SD2.3

SD2.3S



Soft Bone Drills

Code	Description	Price
RSB3	2.8/2.3mmD Step Drill for 3.5mmD Implant	\$60.00
RSB3S	2.8/2.3mmD Short Step Drill for 3.5mmD Implant	\$60.00
RSB4	3.4/2.8mmD Step Drill for 4.3mmD Implant	\$60.00
RSB4S	3.4/2.8mmD Short Step Drill for 4.3mmD Implant	\$60.00
RSB5	4.2/3.3mmD Step Drill for 5.0mmD Implant	\$60.00
RSB5S	4.2/3.3mmD Short Step Drill for 5.0mmD Implant	\$60.00
RSB6	5.1/4.2mmD Step Drill for 6.0mmD Implant	\$60.00
RSB6S	5.1/4.2mmD Short Step Drill for 6.0mmD Implant	\$60.00

RSB3

RSB3S

RSB4

RSB4S



RSB5

RSB5S

RSB6

RSB6S



Dense Bone Drills

Code	Description	Price
RDB3	3.2mmD Dense Bone Drill for 3.5mmD Implant	\$60.00
RDB4	4.0mmD Dense Bone Drill for 4.3mmD Implant	\$60.00
RDB5	4.7mmD Dense Bone Drill for 5.0mmD Implant	\$60.00
RDB6	5.7mmD Dense Bone Drill for 6.0mmD Implant	\$60.00

RDB3

RDB4

RDB5

RDB6



RePlant Insertion Tools

Square connection for use with Ratchet or Torque Wrench

Code	Description	Price
RTT3	Stainless Steel Tri-Lobe Insertion Tool for 3.5mmD Implant	\$60.00
RTT4	Stainless Steel Tri-Lobe Insertion Tool for 4.3mmD Implant	\$60.00
RTT5	Stainless Steel Tri-Lobe Insertion Tool for 5.0mmD Implant	\$60.00
RTT6	Stainless Steel Tri-Lobe Insertion Tool for 6.0mmD Implant	\$60.00

RTT3

RTT4

RTT5

RTT6



RePlant Insertion Drills

Latch head connection for use with Surgical Motor

Code	Description	Price
RTD3	Stainless Steel Tri-Lobe Insertion Drill for 3.5mmD Implant	\$60.00
RTD4	Stainless Steel Tri-Lobe Insertion Drill for 4.3mmD Implant	\$60.00
RTD5	Stainless Steel Tri-Lobe Insertion Drill for 5.0mmD Implant	\$60.00
RTD6	Stainless Steel Tri-Lobe Insertion Drill for 6.0mmD Implant	\$60.00

RTD3

RTD4

RTD5

RTD6





RePlant Drill Sequence

3.5mmD RePlant Straight / Tapered Implant



SD2.3
2.3/2.0mmD

Soft Bone Final



RSB3
2.8/2.3mmD

Dense Bone Final



RDB3
3.2mmD

4.3mmD Straight / Tapered Implant



SD2.3
2.3/2.0mmD



RSB3
2.8/2.3mmD



RSB4
3.4/2.8mmD



RDB4
4.0mmD

5.0mmD Straight / Tapered Implant



SD2.3
2.3/2.0mmD



RSB3
2.8/2.3mmD



RSB4
3.4/2.8mmD



RSB5
4.2/3.3mmD



RDB5
4.7mmD

6.0mmD Straight / Tapered Implant



SD2.3
2.3/2.0mmD



RSB3
2.8/2.3mmD



RSB4
3.4/2.8mmD



RSB5
4.2/3.3mmD




RSB6
5.1/4.2mmD



RDB6
5.7mmD



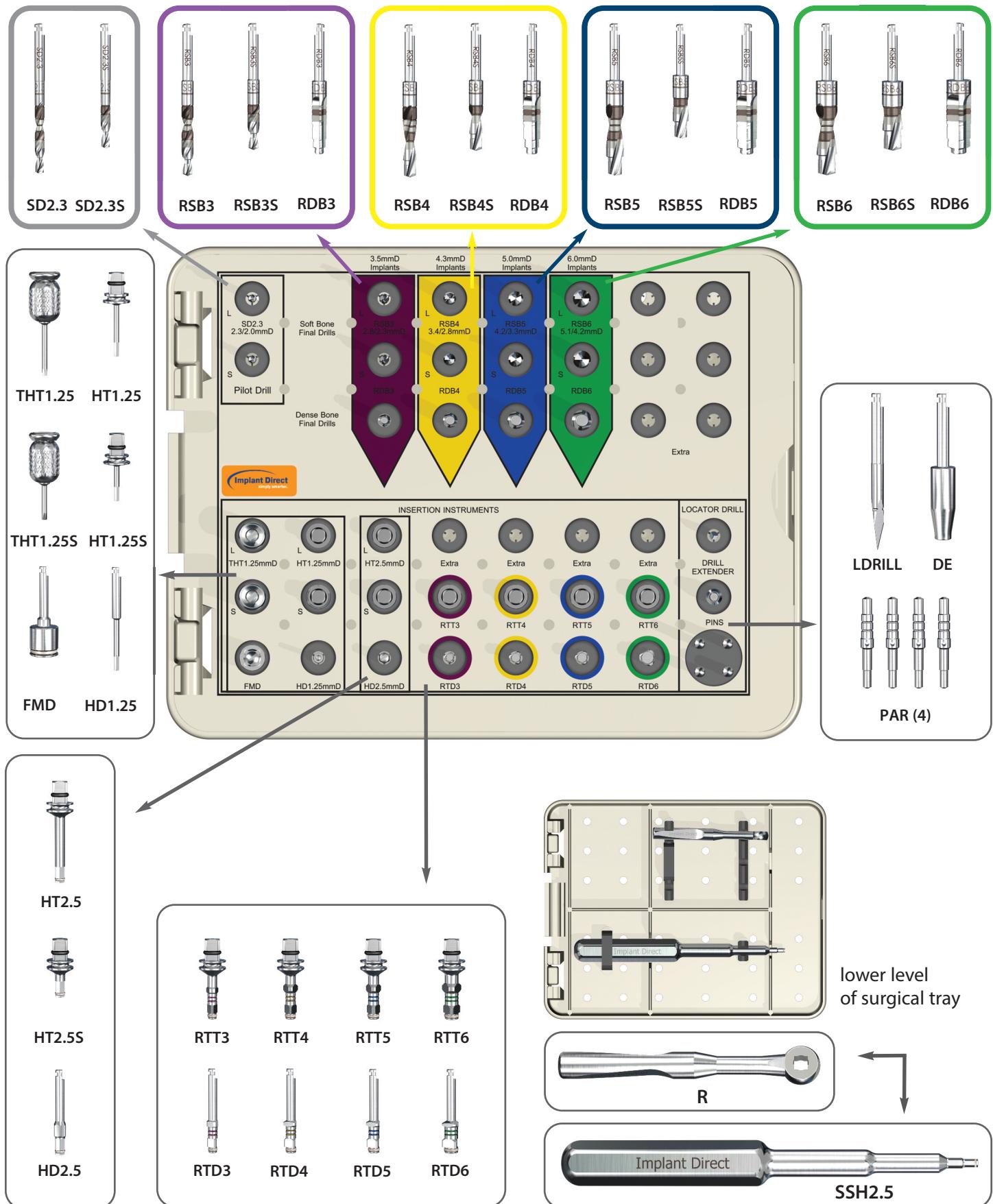Surgical Tray for RePlant® & NobelReplace™ Tapered Groovy Implants






Surgical Trays - RePlant Complete, Basic & Upgrade Kit



Lid



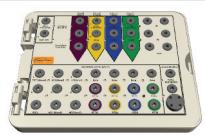Inside

Empty
Tray

Surgical Tray Contents			Complete Kit	Basic Kit	Upgrade Kit
Code	Description	Price	\$2,325.00	\$1,435.00	\$890.00
RST	Empty Surgical Tray	\$250.00	X	X	
SD2.3	2.3/2.0mmD Step Drill-Long	\$60.00	X	X	
SD2.3S	2.3/2.0mmD Step Drill-Short	\$60.00	X		X
RSB3	2.8/2.3mmD Step Drill for 3.5mmD Implant	\$60.00	X	X	
RSB3S	2.8/2.3mmD Short Step Drill for 3.5mmD Implant	\$60.00	X		X
RSB4	3.4/2.8mmD Step Drill for 4.3mmD Implant	\$60.00	X	X	
RSB4S	3.4/2.8mmD Short Step Drill for 4.3mmD Implant	\$60.00	X		X
RSB5	4.2/3.3mmD Step Drill for 5.0mmD Implant	\$60.00	X	X	
RSB5S	4.2/3.3mmD Short Step Drill for 5.0mmD Implant	\$60.00	X		X
RSB6	5.1/4.2mmD Step Drill for 6.0mmD Implant	\$60.00	X	X	
RSB6S	5.1/4.2mmD Short Step Drill for 6.0mmD Implant	\$60.00	X		X
RDB3	3.2mmD Dense Bone Drill for 3.5mmD Implant	\$60.00	X	X	
RDB4	4.0mmD Dense Bone Drill for 4.3mmD Implant	\$60.00	X	X	
RDB5	4.7mmD Dense Bone Drill for 5.0mmD Implant	\$60.00	X	X	
RDB6	5.7mmD Dense Bone Drill for 6.0mmD Implant	\$60.00	X	X	
THT1.25	1.25mmD Tapered Hex Tool-Long	\$40.00	X	X	
THT1.25S	1.25mmD Tapered Hex Tool-Short	\$40.00	X		X
HT1.25	1.25mmD Straight Hex Tool-Long	\$40.00	X	X	
HT1.25S	1.25mmD Straight Hex Tool-Short	\$40.00	X		X
HT2.5	2.5mmD Straight Hex Tool-Long	\$40.00	X		X
HT2.5S	2.5mmD Straight Hex Tool-Short	\$40.00	X	X	
HD1.25	1.25mmD Hex Drill	\$40.00	X	X	
HD2.5	2.5mmD Hex Drill	\$40.00	X		X
RTT3	Tri-Lobe Insertion Tool for 3.5mmD Implant	\$60.00	X	X	
RTT4	Tri-Lobe Insertion Tool for 4.3mmD Implant	\$60.00	X	X	
RTT5	Tri-Lobe Insertion Tool for 5.0mmD Implant	\$60.00	X	X	
RTT6	Tri-Lobe Insertion Tool for 6.0mmD Implant	\$60.00	X	X	
RTD3	Tri-Lobe Insertion Drill for 3.5mmD Implant	\$60.00	X		X
RTD4	Tri-Lobe Insertion Drill for 4.3mmD Implant	\$60.00	X		X
RTD5	Tri-Lobe Insertion Drill for 5.0mmD Implant	\$60.00	X		X
RTD6	Tri-Lobe Insertion Drill for 6.0mmD Implant	\$60.00	X		X
FMD	Fixture Mount Drill	\$80.00	X	X	
LDRILL	Locator Drill	\$30.00	X		X
DE	Drill Extender	\$40.00	X	X	
R	Ratchet	\$125.00	X	X	
PAR (x4)	Paralleling Tool 2.8/2.3	\$60.00	X		X
SSH2.5	2.5mmD Stainless Steel Handle	\$100.00	X		X

Designed For Success, Packaged For Value



Five Year Warranty²

98% Success in
Sinus Elevation Cases³

99% Success in
Immediate Load
Maxillary Cases⁴

Customer
Service &
Technical
Support



Reality Check
Nobel Customers
Save \$518 with ReActive

RePlant® Implant System Advantages:

Nobel Biocare-Compatible Internal Tri-Lobe

Prosthetic compatibility with Nobel Biocare's NobelReplace™ Implant

Surgical Compatibility between RePlant & NobelReplace™

No need to change surgical protocol or tools

Three Implant Designs & Packaging Options

Allows for selection based on price, packaging or thread design

RePlant: \$175 includes cover screw, temporary abutment/transfer

RePlus: \$200 includes cover screw & "snappy" abutment/transfer

ReActive: \$200 includes cover screw & preparable abutment/transfer

Micro-Threads

Reduce crestal stress for improved initial stability

Coronal Bevel on RePlus & ReActive

Allows for use of platform switching technology

Titanium Alloy

Offers superior strength to pure titanium

RePlant®

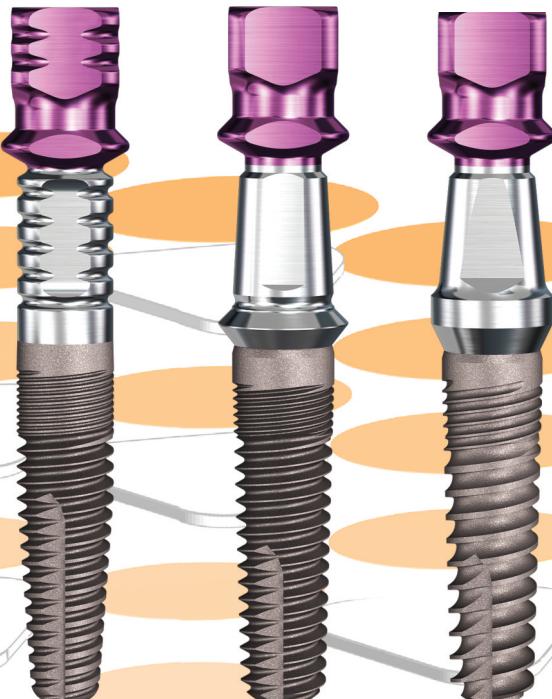
Body Dimensions
Match NobelReplace™

RePlus®

Body Taper
for Bone Expansion

ReActive™

Buttress Threads
for Increased Surface



Laboratory Abutment	Straight Snap-On	Straight Contoured	15° Angled Contoured	Gold/ Zirconia/Ti Plastic Abutment	Plastic Temporary Abutment	Ball GPS™ Attachment	Angled GPS™ Attachment	Multiple-Unit w/Cap & Transfer Straight	Angled
Our price	\$85	\$100	\$100	\$120	\$120	\$40	\$112	\$100	\$120
Nobel Biocare™	\$194	\$192	\$175	\$232	\$223	\$295	\$75	\$225	\$156

Price comparisons based upon US list prices as of January 2013. All trademarks are property of their respective companies.

¹US list price for NobelActive with cover screw, impression coping & abutment

²Terms and conditions apply.

³Implant Dent 2011; 20: 13-19 (n190); ⁴PSN-EAO09PSI-346 (n373) Reports available at www.implantdirect.com

www.implantdirect.com | 888-649-6425



10K-04-14-0044

Printed in the U.S.A. Copyright © 2014 Implant Direct Sybron International. All Rights Reserved